



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

TS4
7-1-2016

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000712 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6-24-16 Cal 1 0.098 AC Cal 2 0.097 AC

Between the dates of 5-26-16 and 6-24-16 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on .
The instrument was returned to service on .

Amy L. Autret

Amy L. Autret


6-29-16

ADPS QAS Permit # 33718

SUBSCRIBED AND SWORN to before me this 29th, day of June, 2016.

Elma Escarcega

Notary Public

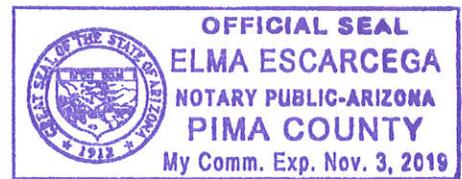

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD QUALITY ASSURANCE PROCEDURE

1 of 8
TSY
7-1-2016

QA SPECIALIST Amy Autret 53546 AGENCY Tucson Police Department

DATE 6-24-2016 TIME 1427

INTOXILYZER SERIAL # 80-000712 LOCATION Crime Laboratory

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS


- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 097 AC

Instrument is operating properly and accurately. YES NO

COMMENTS - Enter key next to apostrophe key. Apostrophe overwrote last name. I wrote in my last name.
- performed additional cal check during SQAP.

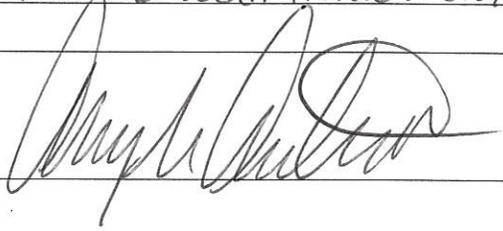SIGNATURE 

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000712
Core Version: 8105.48
06/24/2016

14:44:03

2 of 8
ADD
TS4

Standard Lot#: 06616100A1
Last Changed By: A. AUTRET #53546

QAS: A. 'Autret' #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:45:08
Diagnostic Test	Pass	14:45:42
Air Blank	0.000	14:46:09
0.100 Cal Check	0.097	14:46:30
Air Blank	0.000	14:46:59
Subject Test	0.000	14:47:36
Air Blank	0.000	14:48:03
Wait	ABT*	
Air Blank	0.000	14:50:22

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000712
Core Version: 8105.48
06/24/2016

14:50:39

3068
JL
TJ

Standard Lot#: 06616100A1
Last Changed By: A. AUTRET #53546

QAS: A. 'Autret' #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:51:37
Diagnostic Test	Pass	14:52:11
Air Blank	0.000	14:52:39
0.100 Cal Check	RFI*	14:52:52
Air Blank	IPS**	14:53:00

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000712
Core Version: 8105.48
06/24/2016

14:53:07

4068
All
BY

Standard Lot#: 06616100A1
Last Changed By: A. AUTRET #53546

QAS: A. 'Autret' #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:54:04
Diagnostic Test	Pass	14:54:38
Air Blank	0.000	14:55:06
0.100 Cal Check	0.098	14:55:27
Air Blank	0.000	14:55:56
Subject Test	0.072*	14:59:04
Air Blank	0.000	14:59:36
Air Blank	0.000	15:00:03
Subject Test	XXX**	15:00:29
Air Blank	0.000	15:01:03

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000712
Core Version: 8105.48
06/24/2016 15:01:21

5 of 8
AA
BY

Standard Lot#: 06616100A1
Last Changed By: A. AUTRET #53546

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	15:01:52
0.100 Cal Check	0.097	15:02:13
Air Blank	0.000	15:02:42

Preventative Maintenance
Performed -- Timer Reset

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000712
Core Version: 8105.48
06/24/2016 15:05:48

63F8
JLD
TSG

Standard Lot#: 06616100A1
Last Changed By: A. AUTRET #53546

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	15:06:20
0.100 Cal Check	0.097	15:06:41
Air Blank	0.000	15:07:11

Preventative Maintenance NOT
Performed -- Timer NOT Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

7 of 8
JL
TSY
7-1-2016

QA Specialist Amy Autret 53546 AGENCY Tucson Police Department
DATE 6-24-2016 TIME 1427
Intoxilyzer Serial # 80-000712 Location Crime Laboratory

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 098 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE



DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON PD
Core Version: 8105.48
Serial Number: 80-000712
06/24/2016 14:27:52

Standard Lot#: 06616100A1
Last Changed By: A. AUTRET #53546

QAS: A. ' *Autret* #53546
TUCSON PD

g of 8
dd
TS

Test	g/210L	Time
Air Blank	0.000	14:28:24
0.100 Cal Check	0.098	14:28:45
Air Blank	0.000	14:29:15
End of test		

- ① External printer did not print cal check. Reprinted cal check from database.
- ② Keyboard entry - enter key right next to apostrophe key. Apostrophe overwrote last name. I wrote in my last name where apostrophe overwrote.

Preventive Maintenance ~~NOT~~ *dd*
Performed -- Timer ~~NOT~~ Reset
dd