



CERTIFICATION by CUSTODIAN of RECORDS
 Tucson Police Department Crime Laboratory
 1306 W. Miracle Mile, Tucson, AZ 85705

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-000712 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4
 9.22.2016

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 9/15/2016 Cal 1 0.096 AC Cal 2 0.093 AC

Between the dates of 8/18/2016 and 9/15/2016 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on _____.
 The instrument was returned to service on _____.

William L Loomis

William L. Loomis


9/19/2016

ADPS QAS Permit # 33720

SUBSCRIBED AND SWORN to before me this 19th, day of September, 2016.

Edna Escarcega

Notary Public

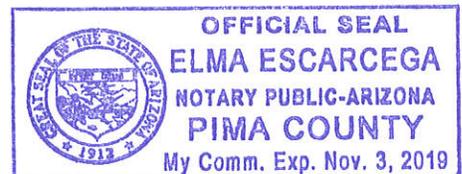

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7
WLL
TS4 9.22.2016

QA SPECIALIST W. LOOMIS 53545 AGENCY Tucson Police Department

DATE 9/15/2016 TIME 13:40

INTOXILYZER SERIAL # 80-000712 LOCATION CRIME LAB

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 093 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Will Loomis

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000712
Core Version: 8105.48
09/15/2016 13:42:03

Standard Lot#: 06616100A1
Last Changed By: A. AUTRET #53546
QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:43:01
Diagnostic Test	Pass	13:43:35
Air Blank	0.000	13:44:02
0.100 Cal Check	0.096	13:44:23
Air Blank	0.000	13:44:52
Subject Test	0.000	13:45:07
Air Blank	0.000	13:45:35
Wait	ABT*	
Air Blank	0.000	13:46:10

*Sequence Aborted

2 of 7
WLL
By

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000712
Core Version: 8105.48
09/15/2016

13:47:43

Standard Lot#: 06616100A1
Last Changed By: A. AUTRET #53546

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:49:02
Diagnostic Test	Pass	13:49:37
Air Blank	0.000	13:50:04
0.100 Cal Check	RFI*	13:50:18
Air Blank	IPS**	13:50:26

*RFI Detect

**Improper Sample

3 of 7
6 WLL
TSy

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000712
Core Version: 8105.48
09/15/2016

13:50:31

Standard Lot#: 06616100A1
Last Changed By: A. AUTRET #53546

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:51:23
Diagnostic Test	Pass	13:51:57
Air Blank	0.000	13:52:25
0.100 Cal Check	0.096	13:52:45
Air Blank	0.000	13:53:15
Subject Test	0.241*	13:56:22
Air Blank	0.000	13:56:56
Air Blank	0.000	13:57:23
Subject Test	XXX**	13:57:40
Air Blank	0.000	13:58:15

*Deficient Sample

**Invalid Sample

Begin new deprivation period

4 of 7
6 ml
TS4

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000712
Core Version: 8105.48
09/15/2016 13:58:32

Standard Lot#: 06616100A1
Last Changed By: A. AUTRET #53546

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	13:59:02
0.100 Cal Check	0.093	13:59:23
Air Blank	0.000	13:59:52

Preventative Maintenance
Performed -- Timer Reset

5 of 7
b wlc

TS4

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD CALIBRATION CHECK PROCEDURE

6 of 7
WLL
TSY 9.22.2016

QA Specialist W. LOOMIS 53545 AGENCY Tucson Police Department
DATE 9/15/2016 TIME 13:40
Intoxilyzer Serial # 80-000712 Location CRIME LAB

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 096 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Will Loomis

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-000712
Core Version: 8105.48
09/15/2016 13:40:29

Standard Lot#: 06616100A1
Last Changed By: A. AUTRET #53546

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	13:41:01
0.100 Cal Check	0.096	13:41:22
Air Blank	0.000	13:41:51

Preventative Maintenance
Performed -- Timer Reset

*7 of 7
WLL
B4*