



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001108 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program. BM

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 2/15/07 Cal 1 0.101 AC Cal 2 0.100 AC

Between the dates of 1/19/07 and 2/15/07 the instrument was operating properly and accurately.


The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Jennifer L. Coombs
Jennifer L. Coombs

2/16/07

ADPS QAS Permit # 33120

SUBSCRIBED AND SWORN to before me this 16, day of February, 2007.

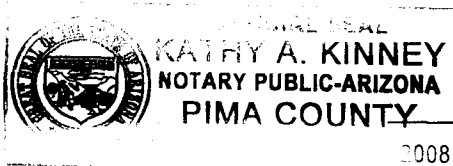

Kathy A. Kinney
Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)1 of 8
file

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE 134

QA SPECIALIST Jennifer Coombs AGENCY Tucson Police DepartmentDATE February 15, 2007 TIME 1402INTOXILYZER SERIAL # 80-001108 LOCATION ODM

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.

(✓) 2. Date check.

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0.000 AC.

(✓) 2. Error recognition logic system functioning.

Not a Successfully Completed Test Sequence printed

(✓) 3. Proper sample recognition system.

Not a Successfully Completed Test Sequence printed ✓✓✓

Deficient sample printed.

(✓) 4. Standard Calibration Check standard 0.100 AC

Result: 0.100 AC

Instrument is operating properly and accurately. YES ☒ NO ☐

COMMENTS

SIGNATURE

Jennifer L. Coombs

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001108
02/15/2007

14:03:52

2 of 8
file
BH

Standard Lot#: 629202I
Last Changed By: J. COOMBS #46160

QAS: J. COOMBS #46160
TUCSON PD CRIME LAB

Operator: J. COOMBS #46160
TUCSON PD CRIME LAB

Subject: SQAP,TEST

DOB: 11/11/1111

Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:04:42
Diagnostic Test	Pass	14:05:10
Air Blank	0.000	14:05:37
0.100 Cal Check	0.101	14:05:56
Air Blank	0.000	14:06:26
Subject Test	0.000	14:06:46
Air Blank	0.000	14:07:13
Five Minute Wait	ABT*	
Air Blank	0.000	14:07:48

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001108
02/15/2007

14:07:53

3 of 8
file
B4

Standard Lot#: 629202I
Last Changed By: J. COOMBS #46160
QAS: J. COOMBS #46160
TUCSON PD CRIME LAB
Operator: J. COOMBS #46160
TUCSON PD CRIME LAB

Subject: SQAP, TEST

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	IPS*	14:08:41
Air Blank	0.000	14:09:09

*Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001108
02/15/2007

14:09:17

418
JL
B4

Standard Lot#: 629202I
Last Changed By: J. COOMBS #46160

QAS: J. COOMBS #46160
TUCSON PD CRIME LAB

Operator: J. COOMBS #46160
TUCSON PD CRIME LAB

Subject: SQAP,TEST

DOB: 11/11/1111

Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:10:34
Diagnostic Test	Pass	14:11:01
Air Blank	0.000	14:11:28
0.100 Cal Check	RFI*	14:11:40
Air Blank	0.000	14:12:11

*RFI Detect

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-001108

02/15/2007

14:12:23

Standard Lot#: 629202I

Last Changed By: J. COOMBS

#46160

QAS: J. COOMBS

#46160

TUCSON PD CRIME LAB

Operator: J. COOMBS

#46160

TUCSON PD CRIME LAB

Subject: SQAP,TEST

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:13:22
Diagnostic Test	Pass	14:13:49
Air Blank	0.000	14:14:17
0.100 Cal Check	0.101	14:14:36
Air Blank	0.000	14:15:06
Subject Test	0.069*	14:18:14
Air Blank	0.000	14:18:46
Air Blank	0.000	14:19:14
Subject Test	XXX**	14:19:33
Air Blank	0.000	14:20:07

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001108

02/15/2007

14:20:24

BH


Standard Lot#: 629202I

Last Changed By: J. COOMBS #46160

QAS: J. COOMBS #46160

TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	14:20:52
Cal Check	0.100	14:21:11
Air Blank	0.000	14:21:41

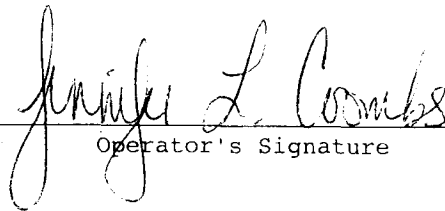
Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Jennifer Coombs 461100 AGENCY Tucson Police Department
DATE February 15, 2007 TIME 1402
Intoxilyzer Serial # 80-001108 Location OSM

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Jennifer L. Coombs

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001108

02/15/2007

14:02:08


Standard Lot#: 629202I

Last Changed By: J. COOMBS #46160

QAS: J. COOMBS #46160

TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	14:02:36
Cal Check	0.101	14:02:54
Air Blank	0.000	14:03:25

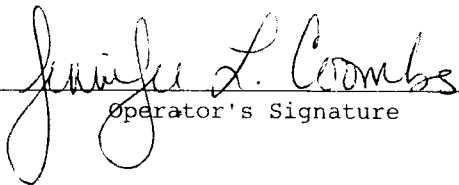

Operator's Signature