



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001108 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4
6-11-07

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6-8-07 Cal 1 0.099 AC Cal 2 0.100 AC

Between the dates of 5-10-07 and 6-8-07
the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.


Clare Conley

Clare L. Conley

6-8-2007

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 8, day of
June, 2007.

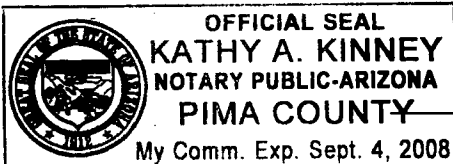

Kathy A Kinney

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

TSy
6-11-07

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Clare Conley 34517 AGENCY Tucson Police Department
DATE 6-8-07 TIME 0655
INTOXILYZER SERIAL # 80-001108 LOCATION ODM

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
(✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES ✓ NO

COMMENTS

SIGNATURE Clare Conley

2 of 8 cu

TS4
6.11.07

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-001108

06/08/2007

06:56:33

Standard Lot#: 629202I

Last Changed By: J. COOMBS

#46160

OAS: C. CONLEY

#34517

TUCSON PD

Operator: C. CONLEY

#34517

TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	06:57:31
Diagnostic Test	Pass	06:57:58
Air Blank	0.000	06:58:25
0.100 Cal Check	0.099	06:58:44
Air Blank	0.000	06:59:14
Subject Test	RFI*	06:59:45
Air Blank	0.000	07:00:13

*RFI Detect

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

3068cu
By
6.11.07

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-001108

06/08/2007

07:00:38

Standard Lot#: 629202I

Last Changed By: J. COOMBS

#46160

QAS: C. CONLEY

#34517

TUCSON PD

Operator: C. CONLEY

#34517

TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:01:32
Diagnostic Test	Pass	07:02:00
Air Blank	0.000	07:02:27
0.100 Cal Check	0.101	07:02:46
Air Blank	0.000	07:03:16
Subject Test	0.000	07:03:34
Air Blank	ABT*	07:03:46
Air Blank	0.000	07:04:14

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

4 of 8 cu
py
6.11.07

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001108
06/08/2007

07:04:56

Standard Lot#: 629202I
Last Changed By: J. COOMBS #46160

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	IPS*	07:06:00
Air Blank	0.000	07:06:28

*Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

5 of 8 cu

TS4
6.11.07

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001108
06/08/2007

07:07:08

Standard Lot#: 629202I
Last Changed By: J. COOMBS #46160
QAS: C. CONLEY #34517
TUCSON PD
Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:08:05
Diagnostic Test	Pass	07:08:32
Air Blank	0.000	07:08:59
0.100 Cal Check	0.099	07:09:18
Air Blank	0.000	07:09:48
Subject Test	0.120*	07:12:56
Air Blank	0.000	07:13:29
Air Blank	0.000	07:13:57
Subject Test	XXX**	07:14:13
Air Blank	0.000	07:14:48

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001108

06/08/2007

07:15:53

Standard Lot#: 629202I

Last Changed By: J. COOMBS #46160

QAS: C. CONLEY #34517

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	07:16:20
Cal Check	0.099	07:16:39
Air Blank	0.000	07:17:09

Clare Conley

Operator's Signature

6 of 8 cu

7/4
6-11-07

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Clare Conley 34517 AGENCY Tucson Police Department
DATE 6-8-07 TIME 0655
Intoxilyzer Serial # 80-001108 Location ODM

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

8 of 8 cu

by
6.11.07

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001108

06/08/2007

06:54:35

Standard Lot#: 629202I

Last Changed By: J. COOMBS #46160

QAS: C. CONLEY #34517

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	06:55:02
Cal Check	0.100	06:55:21
Air Blank	0.000	06:55:51

Clare Conley

Operator's Signature