



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

TS4
8-8-07

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001108 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 8-2-07 Cal 1 0.097 AC Cal 2 0.096 AC

Between the dates of 7-6-07 and 8-2-07 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Clare Conley


Clare L. Conley

8-3-07

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 3, day of

August, 2007.

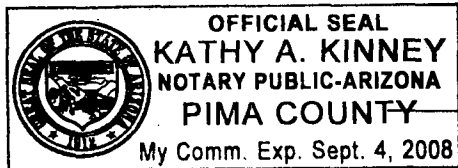

Kathy A Kinney

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000TSY
8.807

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Clare Conley AGENCY Tucson Police DepartmentDATE 8-2-07 TIME 07:08INTOXILYZER SERIAL # 80-001108 LOCATION ODM

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
(✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0.000 AC.
(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(✓) 4. Standard Calibration Check standard 0.100 AC
Result: 0.096 AC

Instrument is operating properly and accurately. YES ☒ NO ☐

COMMENTS Diagnostic test failure (analytical stability
test) occurred due to RFI interference. All
Standard Quality Assurance Procedure checks were normal.

SIGNATURE Clare Conley

269 cu

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-001108

08/02/2007

07:10:58

134

Standard Lot#: 629202I

Last Changed By: J. COOMBS

#46160

QAS: C. CONLEY

#34517

TUCSON PD

Operator: C. CONLEY

#34517

TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:11:50
Diagnostic Test	Pass	07:12:17
Air Blank	0.000	07:12:45
0.100 Cal Check	0.097	07:13:03
Air Blank	0.000	07:13:33
Subject Test	0.000	07:13:49
Air Blank	ABT*	07:13:57
Air Blank	0.000	07:14:25

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

369a

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-001108

08/02/2007

07:15:09

64

Standard Lot#: 629202I

Last Changed By: J. COOMBS

#46160

QAS: C. CONLEY

#34517

TUCSON PD

Operator: C. CONLEY

#34517

TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	07:15:54
Air Blank	0.000	07:16:23

*RFI Detect

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

4896

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-001108

08/02/2007

07:16:28

B4

Standard Lot#: 629202I

Last Changed By: J. COOMBS

#46160

QAS: C. CONLEY

#34517

TUCSON PD

Operator: C. CONLEY

#34517

TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	IPS*	07:16:58
Air Blank	0.000	07:17:27

*Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

5696

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-001108

08/02/2007

07:18:14

B4

Standard Lot#: 629202I

Last Changed By: J. COOMBS

#46160

QAS: C. CONLEY

#34517

TUCSON PD

Operator: C. CONLEY

#34517

TUCSON PD

Subject: TEST,SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:19:09
Diagnostic Test	Fail*	07:19:15
Air Blank	0.000	07:19:43

* Analytical Stability Test

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

60694

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-001108

08/02/2007

07:19:47

BS

Standard Lot#: 629202I

Last Changed By: J. COOMBS

#46160

QAS: C. CONLEY

#34517

TUCSON PD

Operator: C. CONLEY

#34517

TUCSON PD

Subject: TEST,SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:20:44
Diagnostic Test	Pass	07:21:12
Air Blank	0.000	07:21:39
0.100 Cal Check	0.097	07:21:58
Air Blank	0.000	07:22:28
Subject Test	0.021*	07:25:36
Air Blank	0.000	07:26:09
Air Blank	0.000	07:26:36
Subject Test	XXX**	07:26:57
Air Blank	0.000	07:27:33

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 5000

Location: TUCSON POLICE DEPT

Serial Number: 30-001108

08/02/2007

07:27:52

Standard Lot#: 6292021

Last Changed By: J. COOMBS #46160

OAS: C. CONLEY #34517

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	07:28:19
Cal Check	0.096	07:28:38
Air Blank	0.000	07:29:08

Clare Conley

Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

8sf9cu

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

134

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Clare Conley AGENCY Tucson Police Department
DATE 8-2-07 TIME 07:08
Intoxilyzer Serial # 80-001108 Location ODM

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0.100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0.097 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001108

08/02/2007

07:08:56

Standard Lot#: 629202I

Last Changed By: J. COOMBS #46160

QAS: C. CONLEY #34517

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	07:09:23
Cal Check	0.097	07:09:42
Air Blank	0.000	07:10:12

Clare Conley
Operator's Signature