



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001108 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TSy
2.16.2010

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 02/11/10 Cal 1 0.102 AC Cal 2 0.102 AC

Between the dates of 01/13/10 and 02/11/10 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Cheryl K. Straub

Cheryl K. Straub

02/11/10

ADPS QAS Permit # 33408

SUBSCRIBED AND SWORN to before me this 11th day of February, 2010.

A.M. Hernandez

Notary Public



*Includes documents specified in A.A.C. R13-10-104.A.5

197
CS

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD QUALITY ASSURANCE PROCEDURE

TSy
2-16-2010

QA SPECIALIST Cheryl STRAUB AGENCY Tucson Police Department
DATE 02/11/10 TIME 06:05
INTOXILYZER SERIAL # 80-001108 LOCATION TEAM 4

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 102 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Cheryl Straub

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

297
CS

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001108
02/11/2010

06:09:01

124

Standard Lot#: 642845
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	06:10:01
Diagnostic Test	Pass	06:10:28
Air Blank	0.000	06:10:55
0.100 Cal Check	0.103	06:11:14
Air Blank	0.000	06:11:44
Subject Test	0.000	06:12:09
Air Blank	ABT*	06:12:16
Air Blank	0.000	06:12:45

*Sequence Aborted

Not a Successfully Completed Test Sequence

377
CS

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001108
02/11/2010 06:13:07

TS4

Standard Lot#: 642845
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	06:13:45
Air Blank	IPS**	06:13:53

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

497
CS

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001108
02/11/2010

06:14:20

AM

Standard Lot#: 642845
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	06:15:13
Diagnostic Test	Pass	06:15:40
Air Blank	0.000	06:16:07
0.100 Cal Check	0.103	06:16:26
Air Blank	0.000	06:16:56
Subject Test	0.213*	06:20:04
Air Blank	0.000	06:20:37
Air Blank	0.000	06:21:04
Subject Test	XXX**	06:21:29
Air Blank	0.000	06:22:03

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

597
CS


INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001108
02/11/2010 06:22:25

Standard Lot#: 642845
Last Changed By: C. STRAUB #50772

TS4

QAS: C. STRAUB #50772
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	06:22:52
Cal Check	0.102	06:23:11
Air Blank	0.000	06:23:41

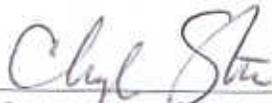

Operator's Signature

EXHIBIT G-2

697
CS

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TSY
2-16-2010

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Cheryl STRAUB AGENCY Tucson Police Department
DATE 02/11/10 TIME 06:05
Intoxilyzer Serial # 80-001108 Location TEAM 4

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 120 AC.
OR
Pour a standard alcohol concentration solution N/A AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 102 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Cheryl Straub

DPS Form Exh G-2 (Rev. 05-01)

787
CS


INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001108
02/11/2010 06:05:28

TS4

Standard Lot#: 642845
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	06:05:56
Cal Check	0.102	06:06:15
Air Blank	0.000	06:06:44

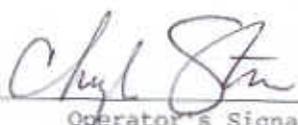

Operator's Signature