



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001108 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

154
5.6.2010

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 05/06/10 Cal 1 0.103 AC Cal 2 0.102 AC

Between the dates of 04/08/10 and 05/06/10 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Cheryl K. Straub

Cheryl K. Straub


05/06/10

ADPS QAS Permit # 33408

SUBSCRIBED AND SWORN to before me this 6th, day of May, 20 10.

A.M. Hernandez

Notary Public

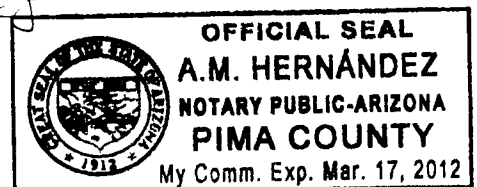

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

177
CS
154
5.6.2010

QA SPECIALIST Cheryl STRAUB AGENCY Tucson Police Department

DATE 05/06/10 TIME 05:58

INTOXILYZER SERIAL # 80-001108 LOCATION TEAM 4

(☒) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (☒) 1. Clock time check.
(☒) 2. Date check.

OPERATIONAL TESTS

- (☒) 1. Alcohol-free subject test result 0. 000 AC.
(☒) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
(☒) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(☒) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 102 AC

Instrument is operating properly and accurately. YES ☒ NO ☐

COMMENTS

SIGNATURE Cheryl Straub

DPS Form Exh G-4 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001108
05/06/2010

06:01:02

277
CS
BY

Standard Lot#: 642845
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 06:02:03 |
| Diagnostic Test | Pass | 06:02:31 |
| Air Blank | 0.000 | 06:02:58 |
| 0.100 Cal Check | 0.103 | 06:03:17 |
| Air Blank | 0.000 | 06:03:46 |
| Subject Test | 0.000 | 06:04:14 |
| Air Blank | ABT* | 06:04:22 |
| Air Blank | 0.000 | 06:04:51 |

*Sequence Aborted

Not a Successfully Completed Test Sequence

377
CS

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001108
05/06/2010

06:04:56

54

Standard Lot#: 642845
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------|--------|----------|
| Air Blank | RFI* | 06:05:43 |
| Air Blank | IPS** | 06:05:54 |

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001108
05/06/2010

06:06:06

487
CS
154

Standard Lot#: 642845
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 06:07:01 |
| Diagnostic Test | Pass | 06:07:28 |
| Air Blank | 0.000 | 06:07:56 |
| 0.100 Cal Check | 0.103 | 06:08:14 |
| Air Blank | 0.000 | 06:08:44 |
| Subject Test | 0.154* | 06:11:52 |
| Air Blank | 0.000 | 06:12:24 |
| Air Blank | 0.000 | 06:12:52 |
| Subject Test | XXX** | 06:13:29 |
| Air Blank | 0.000 | 06:14:03 |

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

577
CS

INTOXILYZER 8000

Location: ARIZONA

Serial Number: 80-001108

05/06/2010

06:14:49

Standard Lot#: 642845

Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772

TUCSON PD

| Test | g/210L | Time |
|-----------|--------|----------|
| Air Blank | 0.000 | 06:15:16 |
| Cal Check | 0.102 | 06:15:35 |
| Air Blank | 0.000 | 06:16:04 |


Operator's Signature

677
CS

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000TS4
5.6.2010

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Cheryl STRAUB AGENCY Tucson Police Department
DATE 05/6/10 TIME 05:58
Intoxilyzer Serial # 80-001108 Location TEAM 4

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution N/A AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 103 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Cheryl Straub

DPS Form Exh G-2 (Rev. 05-01)

777
CS

INTOXILYZER 8000

Location: ARIZONA

Serial Number: 80-001108

05/06/2010

05:58:23

154

Standard Lot#: 642845

Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772

TUCSON PD

| Test | g/210L | Time |
|-----------|--------|----------|
| Air Blank | 0.000 | 05:58:50 |
| Cal Check | 0.103 | 05:59:09 |
| Air Blank | 0.000 | 05:59:38 |



Operator's Signature