



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
270 S. Stone Avenue Tucson, AZ 85701-1917

TJ4  
6-25-2010

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-001108 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6/24/10 Cal 1 0.102 AC Cal 2 0.102 AC

Between the dates of 6/3/10 and 6/24/10 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —  
The instrument was returned to service on —

Sarah Cione

Sarah L. Cione


6/24/10

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 24, day of June, 2010.

AM Hernandez

Notary Public

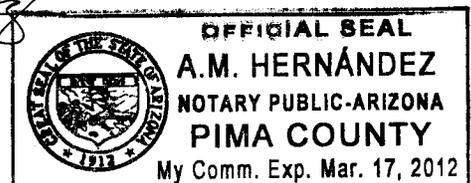

\*Includes documents specified in A.A.C. R13-10-104.A.5

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

TSH  
6.25.2010

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 6/24/10 TIME 1023

INTOXILYZER SERIAL # 80-001108 LOCATION DUI HQ

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 102 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE Sarah Cione

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-001108  
06/24/2010

10:25:27

TS4

---

Standard Lot#: 642845  
Last Changed By: C. STRAUB #50772

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	10:26:13
Diagnostic Test	Pass	10:26:40
Air Blank	0.000	10:27:08
0.100 Cal Check	0.103	10:27:26
Air Blank	0.000	10:27:56
Subject Test	0.000	10:28:15
Air Blank	0.000	10:28:42
Five Minute Wait	ABT*	
Air Blank	0.000	10:29:32

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

3 of 7  
sc

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-001108  
06/24/2010

10:29:38

BS4

---

Standard Lot#: 642845  
Last Changed By: C. STRAUB #50772

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	RFI*	10:30:22
Air Blank	IPS**	10:30:31

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-001108  
06/24/2010

10:30:35

4 of 7 sc

TS4

---

Standard Lot#: 642845  
Last Changed By: C. STRAUB #50772

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	10:31:27
Diagnostic Test	Pass	10:31:54
Air Blank	0.000	10:32:22
0.100 Cal Check	0.103	10:32:41
Air Blank	0.000	10:33:10
Subject Test	0.047*	10:36:18
Air Blank	0.000	10:36:50
Air Blank	0.000	10:37:17
Subject Test	XXX**	10:37:46
Air Blank	0.000	10:38:19

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001108  
06/24/2010 10:38:53

5 of 7 sc

TJH

-----  
Standard Lot#: 642845  
Last Changed By: C. STRAUB #50772  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	10:39:20
Cal Check	0.102	10:39:39
Air Blank	0.000	10:40:08

Sarah Cione  
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

6 of 7 sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

TB4  
6-25-2010

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 6/24/10 TIME 1023  
Intoxilyzer Serial # 80-001108 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 102 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

7 of 7 sc

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001108  
06/24/2010 10:23:22

TS4

-----  
Standard Lot#: 642845  
Last Changed By: C. STRAUB #50772  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	10:23:50
Cal Check	0.102	10:24:09
Air Blank	0.000	10:24:38

Sarah Cione  
Operator's Signature