



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

BY
9.21.2010

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001108 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 9/16/10 Cal 1 0.101 AC Cal 2 0.101 AC

Between the dates of 8/19/10 and 9/16/10 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


9/20/10

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 20th day of September, 2010.

A.M. Hernandez

Notary Public

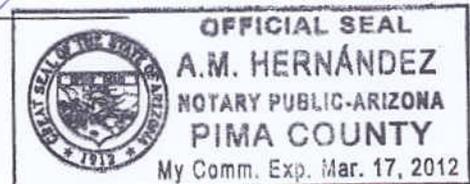

*Includes documents specified in A.A.C. R13-10-104.A.5

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 7 sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TSY
9-21-2010

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 9/16/10 TIME 1126

INTOXILYZER SERIAL # 80-001108 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0.000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0.100 AC
Result: 0.101 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

2 of 7sc

D4

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001108
09/16/2010 11:28:05

Standard Lot#: 642845
Last Changed By: C. STRAUB #50772

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|------------------|--------|----------|
| Air Blank | 0.000 | 11:28:53 |
| Diagnostic Test | Pass | 11:29:20 |
| Air Blank | 0.000 | 11:29:48 |
| 0.100 Cal Check | 0.101 | 11:30:06 |
| Air Blank | 0.000 | 11:30:36 |
| Subject Test | 0.000 | 11:31:05 |
| Air Blank | 0.000 | 11:31:33 |
| Five Minute Wait | ABT* | |
| Air Blank | 0.000 | 11:32:38 |

*Sequence Aborted

Not a Successfully Completed Test Sequence

3 of 7 sc
B4

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001108
09/16/2010 11:33:03

Standard Lot#: 642845
Last Changed By: C. STRAUB #50772

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------|--------|----------|
| Air Blank | RFI* | 11:33:30 |
| Air Blank | IPS** | 11:33:37 |

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001108
09/16/2010

11:35:32

4 of 7 sc

154

Standard Lot#: 642845
Last Changed By: C. STRAUB #50772

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 11:36:41 |
| Diagnostic Test | Pass | 11:37:09 |
| Air Blank | 0.000 | 11:37:36 |
| 0.100 Cal Check | 0.101 | 11:37:55 |
| Air Blank | 0.000 | 11:38:24 |
| Subject Test | 0.159* | 11:41:33 |
| Air Blank | 0.000 | 11:42:04 |
| Air Blank | 0.000 | 11:42:32 |
| Subject Test | XXX** | 11:42:54 |
| Air Blank | 0.000 | 11:43:27 |

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001108
09/16/2010 11:43:46

5 of 7 SC

TS4

Standard Lot#: 642845
Last Changed By: C. STRAUB #50772

QAS: S. CIONE #50616
TUCSON PD

| Test | g/210L | Time |
|-----------|--------|----------|
| Air Blank | 0.000 | 11:44:13 |
| Cal Check | 0.101 | 11:44:32 |
| Air Blank | 0.000 | 11:45:02 |

Sarah Cione
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

6 of 7 SC

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

By
9.21.2010

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 9/16/10 TIME 1126
Intoxilyzer Serial # 80-001108 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001108
09/16/2010

11:26:12

7 of 7 sc

BY

Standard Lot#: 642845
Last Changed By: C. STRAUB #50772

QAS: S. CIONE #50616
TUCSON PD

| Test | g/210L | Time |
|-----------|--------|----------|
| Air Blank | 0.000 | 11:26:39 |
| Cal Check | 0.101 | 11:26:58 |
| Air Blank | 0.000 | 11:27:27 |

Sarah Cione
Operator's Signature