



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

TS4
10-18-2010

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001108 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 10-14-10 Cal 1 0.101 AC Cal 2 0.101 AC

Between the dates of 9-16-10 and 10-14-10 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Clare Conley

Clare L. Conley


10-15-10

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 15, day of October, 2010.

Laura B. Hodgson

Notary Public

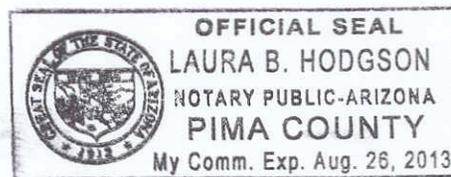

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1587
Cu
TS4
10.18.2010

QA SPECIALIST Clare Conley AGENCY Tucson Police Department

DATE 10-14-10 TIME 0932

INTOXILYZER SERIAL # 80-001108 LOCATION ODWJUH

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 101 AC

Instrument is operating properly and accurately. YES NO

COMMENTS None

SIGNATURE Clare Conley

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001108
10/14/2010

09:35:20

287
CW
154

Standard Lot#: 642845
Last Changed By: C. STRAUB #50772

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:36:07
Diagnostic Test	Pass	09:36:35
Air Blank	0.000	09:37:02
0.100 Cal Check	0.101	09:37:21
Air Blank	0.000	09:37:51
Subject Test	0.000	09:38:27
Air Blank	ABT*	09:38:34
Air Blank	0.000	09:39:03

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001108
10/14/2010

09:39:18

3A7
CM
184

Standard Lot#: 642845
Last Changed By: C. STRAUB #50772

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:40:11
Diagnostic Test	Pass	09:40:39
Air Blank	0.000	09:41:06
0.100 Cal Check	0.100	09:41:25
Air Blank	0.000	09:41:54
Subject Test	RFI*	09:42:06
Air Blank	IPS**	09:42:27

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001108
10/14/2010

09:42:36

4567
Cu
101

Standard Lot#: 642845
Last Changed By: C. STRAUB #50772

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:43:23
Diagnostic Test	Pass	09:43:51
Air Blank	0.000	09:44:18
0.100 Cal Check	0.101	09:44:37
Air Blank	0.000	09:45:06
Subject Test	0.157*	09:48:14
Air Blank	0.000	09:48:46
Air Blank	0.000	09:49:14
Subject Test	XXX**	09:50:01
Air Blank	0.000	09:50:34

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001108
10/14/2010

09:51:44

Standard Lot#: 642845
Last Changed By: C. STRAUB #50772

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:52:11
Cal Check	0.101	09:52:30
Air Blank	0.000	09:52:59

Clare Conley

Operator's Signature

507
W
154

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

b of 7
CW

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TS4
10.18.2010

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Clare Conley AGENCY Tucson Police Department
DATE 10-14-10 TIME 0932
Intoxilyzer Serial # 80-00 1108 Location DDW DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001108
10/14/2010 09:32:34


787
CW

Standard Lot#: 642845
Last Changed By: C. STRAUB #50772

QAS: C. CONLEY #34517
TUCSON PD

154

Test	g/210L	Time
Air Blank	0.000	09:33:01
Cal Check	0.101	09:33:20
Air Blank	0.000	09:33:50

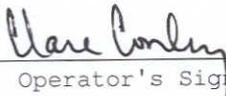

Operator's Signature