



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-001108 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

154  
2-7-2011

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 2-4-11 Cal 1 0.102 AC Cal 2 0.101 AC

Between the dates of 1-6-11 and 2-4-11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.  
The instrument was returned to service on N/A.


Clare Conley

Clare L. Conley

2-4-11

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 4th, day of February, 2011.

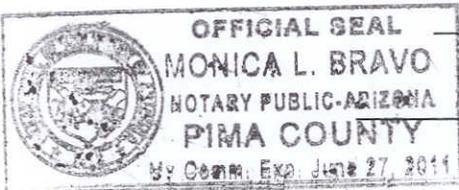

Monica L. Bravo  
Monica L. Bravo  
Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

P1067  
CM

TS4  
2-7-2011

QA SPECIALIST Clare Conley AGENCY Tucson Police Department  
DATE 2-4-11 TIME 0949  
INTOXILYZER SERIAL # 80-001108 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 101 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS None  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Clare Conley  
\_\_\_\_\_

2067  
cu

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-001108  
02/04/2011

09:52:54

SM

---

Standard Lot#: 20910100A1  
Last Changed By: C. CONLEY #34517

---

QAS: C. CONLEY #34517  
TUCSON PD

---

Operator: C. CONLEY #34517  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	09:53:39
Diagnostic Test	Pass	09:54:07
Air Blank	0.000	09:54:34
0.100 Cal Check	0.102	09:54:53
Air Blank	0.000	09:55:22
Subject Test	0.000	09:55:52
Air Blank	ABT*	09:55:59
Air Blank	0.000	09:56:27

\*Sequence Aborted

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

3567  
CW

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-001108  
02/04/2011 09:56:32

BH

---

Standard Lot#: 20910100A1  
Last Changed By: C. CONLEY #34517

---

QAS: C. CONLEY #34517  
TUCSON PD

---

Operator: C. CONLEY #34517  
TUCSON PD

---

Subject: TEST,SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	IPS*	09:57:09
Air Blank	RFI**	09:57:17

\*Improper Sample  
\*\*RFI Detect

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

4067  
cu

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-001108  
02/04/2011 09:57:24

134

---

Standard Lot#: 20910100A1  
Last Changed By: C. CONLEY #34517

---

QAS: C. CONLEY #34517  
TUCSON PD

---

Operator: C. CONLEY #34517  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	09:58:11
Diagnostic Test	Pass	09:58:39
Air Blank	0.000	09:59:06
0.100 Cal Check	0.102	09:59:25
Air Blank	0.000	09:59:54
Subject Test	0.240*	10:03:03
Air Blank	0.000	10:03:35
Air Blank	0.000	10:04:02
Subject Test	XXX**	10:04:28
Air Blank	0.000	10:05:01

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001108  
02/04/2011

10:05:17

-----  
Standard Lot#: 20910100A1  
Last Changed By: C. CONLEY #34517  
-----

QAS: C. CONLEY #34517  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	10:05:45
Cal Check	0.101	10:06:03
Air Blank	0.000	10:06:33

*Clare Conley*

\_\_\_\_\_  
Operator's Signature

5067  
CM

BY

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

b5f7  
cu

TSY  
2.7.2011

QA Specialist Clare Conley AGENCY Tucson Police Department  
DATE 2-4-11 TIME 0949  
Intoxilyzer Serial # 80-001108 Location Crime Lab

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 102 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001108  
02/04/2011 09:49:03

7 of 7  
en

-----  
Standard Lot#: 20910100A1  
Last Changed By: C. CONLEY #34517  
-----

154

-----  
QAS: C. CONLEY #34517  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:49:31
Cal Check	0.102	09:49:49
Air Blank	0.000	09:50:19

*Clare Conley*

\_\_\_\_\_  
Operator's Signature