



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001108 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4
6-27-2011

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6/23/11 Cal 1 0.102 AC Cal 2 0.102 AC

Between the dates of 5/26/11 and 6/23/11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -
The instrument was returned to service on -

Sarah Cione

Sarah L. Cione


6/27/11

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 27, day of June, 2011.

A.M. Hernandez

Notary Public

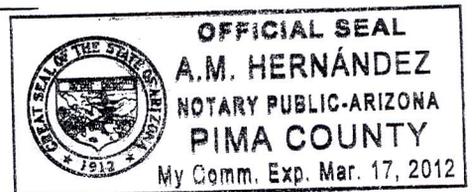

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 7 SC

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TSY
6.27.2011

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 6/23/11 TIME 1352

INTOXILYZER SERIAL # 80-001108 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 102 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001108
06/23/2011

13:54:24

2 of 7sc

134

Standard Lot#: 20910100A1
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:55:13
Diagnostic Test	Pass	13:55:41
Air Blank	0.000	13:56:08
0.100 Cal Check	0.102	13:56:27
Air Blank	0.000	13:56:56
Subject Test	0.000	13:57:15
Air Blank	0.000	13:57:42
Five Minute Wait	ABT*	
Air Blank	0.000	13:58:18

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001108
06/23/2011

13:58:57

3 of 7_{sc}

BSH

Standard Lot#: 20910100A1
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	13:59:29
Air Blank	IPS**	13:59:37

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001108
06/23/2011

13:59:43

4 of 7 sc
154

Standard Lot#: 20910100A1
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:00:32
Diagnostic Test	Pass	14:01:00
Air Blank	0.000	14:01:27
0.100 Cal Check	0.103	14:01:46
Air Blank	0.000	14:02:15
Subject Test	0.038*	14:05:24
Air Blank	0.000	14:05:55
Air Blank	0.000	14:06:22
Subject Test	XXX**	14:06:45
Air Blank	0.000	14:07:18

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001108
06/23/2011 14:07:46

5 of 7 sc

Standard Lot#: 20910100A1
Last Changed By: C. CONLEY #34517

TS4

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:08:14
Cal Check	0.102	14:08:33
Air Blank	0.000	14:09:02

Sarah Cione

Operator's Signature

EXHIBIT G-2

6 of 7 sc

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD CALIBRATION CHECK PROCEDURE

T54
6-27-2011

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 6/23/11 TIME 1352
Intoxilyzer Serial # 80-001108 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 102 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001108
06/23/2011 13:52:51

7 of 7 SC

B4

Standard Lot#: 20910100A1
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	13:53:18
Cal Check	0.102	13:53:37
Air Blank	0.000	13:54:06

Sarah Cione

Operator's Signature