



# CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

769  
9.28.2011

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-001108 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 9/15/11 Cal 1 0.099 AC Cal 2 0.099 AC

Between the dates of 8/18/11 and 9/15/11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


9/22/11

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 22, day of September, 2011.

A.M. Hernandez

Notary Public

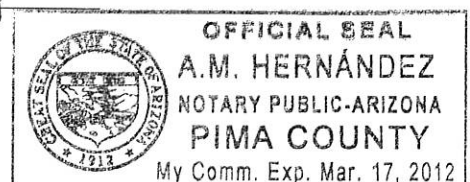

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7<sub>sc</sub>  
By  
9.28.2011

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 9/15/11 TIME 1321

INTOXILYZER SERIAL # 80-001108 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.

(✓) 2. Date check.

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0. 000 AC.

(✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed

(✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.

(✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 099 AC

Instrument is operating properly and accurately. YES ☒ NO ☐

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Sarah Cione

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-001108  
09/15/2011

13:23:49

Standard Lot#: 20910100A1  
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:24:59
Diagnostic Test	Pass	13:25:26
Air Blank	0.000	13:25:53
0.100 Cal Check	0.100	13:26:12
Air Blank	0.000	13:26:42
Subject Test	0.000	13:27:09
Air Blank	0.000	13:27:36
Five Minute Wait	ABT*	
Air Blank	0.000	13:28:37

\*Sequence Aborted

Not a Successfully Completed Test Sequence

3 of 7<sub>sc</sub>

154

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-001108  
09/15/2011

13:28:46

---

Standard Lot#: 20910100A1  
Last Changed By: C. CONLEY #34517

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

---

15 Minute Deprivation Period? Yes

---

---

Test	g/210L	Time
Air Blank	RFI*	13:29:34
Air Blank	IPS**	13:29:43

---

\*RFI Detect

\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

4 of 7<sub>sc</sub>

BH

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-001108  
09/15/2011 13:29:52

Standard Lot#: 20910100A1  
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:30:44
Diagnostic Test	Pass	13:31:12
Air Blank	0.000	13:31:39
0.100 Cal Check	0.101	13:31:58
Air Blank	0.000	13:32:28
Subject Test	0.036*	13:35:36
Air Blank	0.000	13:36:09
Air Blank	0.000	13:36:36
Subject Test	XXX**	13:36:52
Air Blank	0.000	13:37:25

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

5 of 7<sub>sc</sub>

B4

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001108  
09/15/2011 13:37:51

-----  
Standard Lot#: 20910100A1  
Last Changed By: C. CONLEY #34517  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	13:38:18
Cal Check	0.099	13:38:37
Air Blank	0.000	13:39:06

  
\_\_\_\_\_  
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 9/15/11 TIME 1321  
Intoxilyzer Serial # 80-001108 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach  $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: ARIZONA

Serial Number: 80-001108

09/15/2011

13:21:14

-----  
Standard Lot#: 20910100A1

Last Changed By: C. CONLEY #34517  
-----

QAS: S. CIONE #50616

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	13:21:42
Cal Check	0.099	13:22:01
Air Blank	0.000	13:22:30

*Sarah Cione*

\_\_\_\_\_  
Operator's Signature

7 of 7sc  
154