



CERTIFICATION by CUSTODIAN of RECORDS  
Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-001108 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4  
11.7-2011

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 10/14/11 Cal 1 0.101 AC Cal 2 0.101 AC

Between the dates of 9/15/11 and 10/14/11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


10/27/11

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 27, day of October, 2011.

A.M. Hernandez

Notary Public

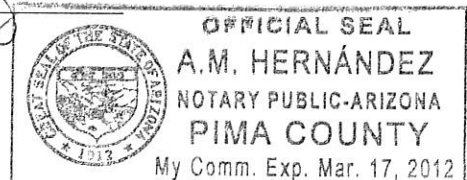

\*Includes documents specified in A.A.C. R13-10-104.A.5

## EXHIBIT G-4

1078<sub>sc</sub>

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

TS4  
11-7-2011QA SPECIALIST Sarah Cione AGENCY Tucson Police DepartmentDATE 10/14/11 TIME 1414INTOXILYZER SERIAL # 80-001108 LOCATION CrimeLab

(✓) 1. Display reads "PUSH BUTTON TO START"

## DIAGNOSTIC TESTS

(✓) 1. Clock time check.

(✓) 2. Date check.

## OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0. 000 AC.(✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed(✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.(✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 101 ACInstrument is operating properly and accurately. YES ✓ NO       

## COMMENTS

---

---

---

---

SIGNATURE Sarah Cione

2 of 8<sub>sc</sub>

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-001108  
10/14/2011

14:18:20

TJ4

Standard Lot#: 20910100A1  
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test             | g/210L | Time     |
|------------------|--------|----------|
| Air Blank        | 0.000  | 14:19:10 |
| Diagnostic Test  | Pass   | 14:19:38 |
| Air Blank        | 0.000  | 14:20:05 |
| 0.100 Cal Check  | 0.100  | 14:20:24 |
| Air Blank        | 0.000  | 14:20:53 |
| Subject Test     | 0.000  | 14:21:15 |
| Air Blank        | 0.000  | 14:21:42 |
| Five Minute Wait | ABT*   |          |
| Air Blank        | 0.000  | 14:22:20 |

\*Sequence Aborted

Not a Successfully Completed Test Sequence

3 of 8 sc

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-001108  
10/14/2011

14:23:23

TSH

Standard Lot#: 20910100A1  
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111

Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test      | g/210L | Time     |
|-----------|--------|----------|
| Air Blank | RFI*   | 14:23:54 |
| Air Blank | IPS**  | 14:24:02 |

\*RFI Detect

\*\*Improper Sample

Not a Successfully Completed Test Sequence

4 of 8<sub>sc</sub>

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-001108  
10/14/2011

14:24:50

TS4

Standard Lot#: 20910100A1  
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

| Test            | g/210L | Time     |
|-----------------|--------|----------|
| Air Blank       | 0.000  | 14:25:37 |
| Diagnostic Test | Fail*  | 14:25:41 |
| Air Blank       | 0.000  | 14:26:09 |

\* Voltage/Current Test

Not a Successfully Completed Test Sequence

5 of 8<sub>sc</sub>

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-001108  
10/14/2011

14:26:14

B4

Standard Lot#: 20910100A1  
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

| Test            | g/210L | Time     |
|-----------------|--------|----------|
| Air Blank       | 0.000  | 14:27:03 |
| Diagnostic Test | Pass   | 14:27:30 |
| Air Blank       | 0.000  | 14:27:58 |
| 0.100 Cal Check | 0.102  | 14:28:16 |
| Air Blank       | 0.000  | 14:28:46 |
| Subject Test    | 0.036* | 14:31:54 |
| Air Blank       | 0.000  | 14:32:25 |
| Air Blank       | 0.000  | 14:32:53 |
| Subject Test    | XXX**  | 14:33:15 |
| Air Blank       | 0.000  | 14:33:48 |

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

6 of 8<sub>sc</sub>

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001108  
10/14/2011 14:34:05

TS4

-----  
Standard Lot#: 20910100A1  
Last Changed By: C. CONLEY #34517  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

| Test      | g/210L | Time     |
|-----------|--------|----------|
| Air Blank | 0.000  | 14:34:32 |
| Cal Check | 0.101  | 14:34:51 |
| Air Blank | 0.000  | 14:35:20 |

Sarah Cione

Operator's Signature

## EXHIBIT G-2

7 of 8 sc

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

TS4  
11.7.2011

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 10/14/11 TIME 1414  
Intoxilyzer Serial # 80-001108 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach  $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)



INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001108  
10/14/2011

14:14:33

TS4

-----  
Standard Lot#: 20910100A1  
Last Changed By: C. CONLEY #34517  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

| Test      | g/210L | Time     |
|-----------|--------|----------|
| Air Blank | 0.000  | 14:15:01 |
| Cal Check | 0.101  | 14:15:19 |
| Air Blank | 0.000  | 14:15:49 |

Sarah Cione  
Operator's Signature