



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

TS4  
6-12-2012

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-001108 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 5/24/12 Cal 1 0.101 AC Cal 2 0.100 AC

Between the dates of 4/26/12 and 5/24/12 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


5/24/12

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 24th, day of May, 2012.

Elma Escarcega

Notary Public

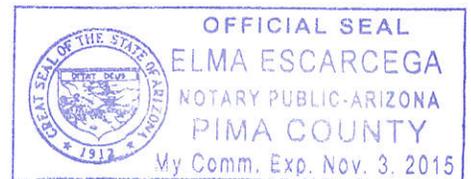

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

1 of 7 sc  
TSh  
6-12-2012

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 5/24/12 TIME 927

INTOXILYZER SERIAL # 80-001108 LOCATION DUI HQ

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0.000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0.100 AC  
Result: 0.100 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

2 of 7  
5c  
PH

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-001108  
Core Version: 8105.48  
05/24/2012 09:29:35

---

Standard Lot#: 20910100A1  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

| Test            | g/210L | Time     |
|-----------------|--------|----------|
| Air Blank       | 0.000  | 09:30:23 |
| Diagnostic Test | Pass   | 09:30:58 |
| Air Blank       | 0.000  | 09:31:26 |
| 0.100 Cal Check | 0.101  | 09:31:47 |
| Air Blank       | 0.000  | 09:32:16 |
| Subject Test    | 0.000  | 09:32:31 |
| Air Blank       | 0.000  | 09:32:58 |
| Wait            | ABT*   |          |
| Air Blank       | 0.000  | 09:34:20 |

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

3 of 7 sc  
154

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-001108  
Core Version: 8105.48  
05/24/2012 09:34:54

---

Standard Lot#: 20910100A1  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

| Test      | g/210L | Time     |
|-----------|--------|----------|
| Air Blank | RFI*   | 09:35:20 |
| Air Blank | IPS**  | 09:35:28 |

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

4 of 7 sc  
T24

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-001108  
Core Version: 8105.48  
05/24/2012 09:36:57

---

Standard Lot#: 20910100A1  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

| Test            | g/210L | Time     |
|-----------------|--------|----------|
| Air Blank       | 0.000  | 09:37:39 |
| Diagnostic Test | Pass   | 09:38:14 |
| Air Blank       | 0.000  | 09:38:41 |
| 0.100 Cal Check | 0.101  | 09:39:02 |
| Air Blank       | 0.000  | 09:39:31 |
| Subject Test    | 0.053* | 09:42:39 |
| Air Blank       | 0.000  | 09:43:11 |
| Air Blank       | 0.000  | 09:43:38 |
| Subject Test    | XXX**  | 09:43:57 |
| Air Blank       | 0.000  | 09:44:30 |

\*Deficient Sample  
\*\*Invalid Sample  
Begin new deprivation period

**Not a Successfully Completed Test Sequence**

5 of 7 sc  
T54

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-001108  
Core Version: 8105.48  
05/24/2012 09:44:51

-----  
Standard Lot#: 20910100A1  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

| Test            | g/210L | Time     |
|-----------------|--------|----------|
| Air Blank       | 0.000  | 09:45:19 |
| 0.100 Cal Check | 0.100  | 09:45:40 |
| Air Blank       | 0.000  | 09:46:10 |

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

6 of 7<sub>sc</sub>

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

TS4  
6.12.2012

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 5/24/12 TIME 927  
Intoxilyzer Serial # 80-001108 Location DU.I HQ

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 101 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-001108  
Core Version: 8105.48  
05/24/2012 09:27:19

7 of 7 sc  
134

-----  
Standard Lot#: 20910100A1  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

| Test            | g/210L | Time     |
|-----------------|--------|----------|
| Air Blank       | 0.000  | 09:27:49 |
| 0.100 Cal Check | 0.101  | 09:28:10 |
| Air Blank       | 0.000  | 09:28:39 |

Preventative Maintenance  
Performed -- Timer Reset