



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

TS4  
9.4.2014

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-001108 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 8/15/14 Cal 1 0.097 AC Cal 2 0.096 AC

Between the dates of 7/18/14 and 8/15/14 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione

8/19/14


ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 19th, day of August, 2014.

Elma Escarcega

Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

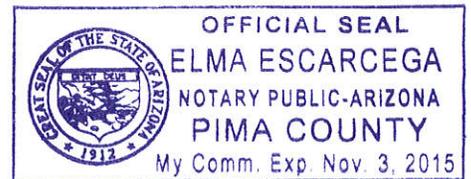

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7  
sc  
TS4  
9.4.2014

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 8/15/14 TIME 1620

INTOXILYZER SERIAL # 80-001108 LOCATION Crime Lab

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 096 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Sarah Cione

2 of 7  
Σ

T34

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-001108  
Core Version: 8105.48  
08/15/2014 16:22:26

-----  
Standard Lot#: 25913100A2  
Last Changed By: S. CIONE #50616  
-----

DAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Subject: SOAP,SOAP  
DOB: 11/11/1980  
Sex: F Weight: 111  
-----

15 Minute Deprivation Period? Yes  
-----

Test	g/210L	Time
Air Blank	0.000	16:23:23
Diagnostic Test	Pass	16:24:01
Air Blank	0.000	16:24:28
0.100 Cal Check	0.097	16:24:49
Air Blank	0.000	16:25:19
Subject Test	0.000	16:25:34
Air Blank	0.000	16:26:01
Wait	ABT*	
Air Blank	0.000	16:26:45

\*Sequence Aborted

Not a Successfully  
Completed  
Test Sequence

3 of 7 sc  
TS4

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-001108  
Core Version: 8105.48  
08/15/2014 16:27:11

Standard Lot#: 25913100A2  
Last Changed By: S. CIONE #50616

OAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SOAP, SOAP  
DOB: 11/11/1980  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:27:53
Diagnostic Test	Pass	16:28:30
Air Blank	0.000	16:28:58
0.100 Cal Check	0.097	16:29:19
Air Blank	0.000	16:29:48
Subject Test	RFI*	16:30:07
Air Blank	IPS**	16:30:15

\*RFI Detect  
\*\*Improper Sample

Not a Successfully  
Completed  
Test Sequence

4 of 7  
r sc

T34

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-001108  
Core Version: 8105.48  
08/15/2014 16:30:29

Standard Lot#: 25913100A2  
Last Changed By: S. CIONE #50616

OAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SOAP, SOAP  
DOB: 11/11/1980  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:31:10
Diagnostic Test	Pass	16:31:47
Air Blank	0.000	16:32:14
0.100 Cal Check	0.097	16:32:35
Air Blank	0.000	16:33:05
Subject Test	0.260*	16:36:12
Air Blank	0.000	16:36:45
Air Blank	0.000	16:37:12
Subject Test	XXX**	16:37:37
Air Blank	0.000	16:38:10

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

Not a Successfully  
Completed  
Test Sequence

5 of 7<sub>sc</sub>

734

---

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-001108  
Core Version: 8105.48  
08/15/2014 16:38:37

---

Standard Lot#: 25913100A2  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Test	g/210L	Time
Air Blank	0.000	16:39:05
0.100 Cal Check	0.096	16:39:26
Air Blank	0.000	16:39:56

---

Preventative Maintenance  
Performed -- Timer Reset

---

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7  
5c  
TS4  
9.4.2014

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 8/15/14 TIME 1620  
Intoxilyzer Serial # 80-001108 Location Crime Lab

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 097 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

7 of 7  
5c  
TS4

---

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-001108  
Core Version: 8105.48  
08/15/2014 16:20:03

---

Standard Lot#: 25913100A2  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Test	g/210L	Time
Air Blank	0.000	16:20:31
0.100 Cal Check	0.097	16:20:52
Air Blank	0.000	16:21:22

---

Preventative Maintenance  
Performed -- Timer Reset

---