



CERTIFICATION by CUSTODIAN of RECORDS  
 Tucson Police Department Crime Laboratory  
 1306 W. Miracle Mile, Tucson, AZ 85705

TS4  
 3.11.2015

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-001108 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 3-5-15 Cal 1 0.097 AC Cal 2 0.098 AC

Between the dates of 2-5-15 and 3-5-15 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.  
 The instrument was returned to service on N/A.

*Amy L. Autret*  
 Amy L. Autret


3-9-15

ADPS QAS Permit # 33718

SUBSCRIBED AND SWORN to before me this 9<sup>th</sup>, day of March, 2015.

*Elma Escarcega*

Notary Public

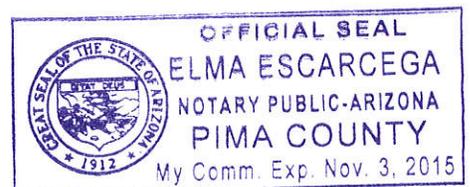

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

1 of 8  
JLH  
TBJ  
3.11.2015

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Amy Autret 53546 AGENCY Tucson Police Department

DATE 3/5/15 TIME 1403

INTOXILYZER SERIAL # 80-001108 LOCATION Lab

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS


- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 098 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE

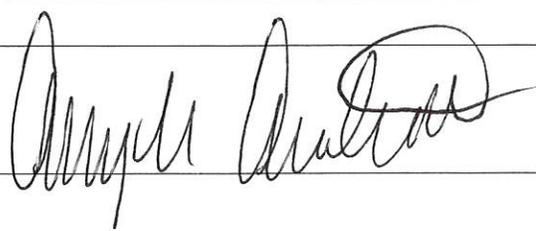

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

2 of 8  
MA  
Tdy

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-001108  
Core Version: 8105.48  
03/05/2015 14:10:28

Standard Lot#: 2591310042  
Last Changed By: G. CLONE #50616

QAS: A. AUTRET #53546  
TUCSON PD

Operator: A. AUTRET #53546  
TUCSON PD

Subject: SOAP, SOAP  
DOB: 11/11/1987  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:11:25
Diagnostic Test	Pass	14:12:03
Air Blank	0.000	14:12:30
0.100 Cal Check	0.098	14:12:51
Air Blank	0.000	14:13:20
Subject Test	0.000	14:13:37
Air Blank	0.000	14:14:04
Wait	ABT*	
Air Blank	0.000	14:14:40

\*Sequence Aborted

Not a Successfully  
Completed  
Test Sequence

3 of 8  
AB  
TS4

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-001108  
Core Version: 8105.48  
03/05/2015 14:15:41

Standard Lot#: 25913100A2  
Last Changed By: S. CLONE #53616

OAS: A. AUTRET #53546  
TUCSON PD

Operator: A. AUTRET #53546  
TUCSON PD

Subject: SOAP, SOAP  
DOB: 11/11/1987  
Sex: F Weight: 111

15 Minute Depriuation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:16:50
Diagnostic Test	Pass	14:17:28
Air Blank	0.000	14:17:55
0.100 Cal Check	RFI*	14:18:07
Air Blank	IPS**	14:18:15

\*RFI Detect  
\*\*Improper Sample

Not a Successfully  
Completed  
Test Sequence

4 of 8  
AD  
TS4

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-001108  
Core Version: 8105.48  
03/05/2015 14:19:00

Standard Length: 25913100A2  
Last Changed By: S. CIONE #50616

OAS: A. AUTRET #53546  
TUCSON PD

Operator: A. AUTRET #53546  
TUCSON PD

Subject: SAP, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:20:08
Diagnostic Test	Pass	14:20:46
Air Blank	0.000	14:21:13
0.100 Cal Check	0.099	14:21:34
Air Blank	0.000	14:22:03
Subject Test	0.089*	14:25:11
Air Blank	0.000	14:25:42
Air Blank	0.000	14:26:10
Subject Test	0.167**	14:26:34
Air Blank	0.000	14:27:07
Wait		
Air Blank	0.000	14:32:02
Subject Test	0.000**	14:32:19
Air Blank	0.000	14:32:47

\*Deficient Sample  
\*\*No .02 agreement

Run out of breath  
for Invalid sample.  
waited 5 minutes,  
0.000 result.

Not a Successfully  
Completed  
Test Sequence

5 of 8  
AS  
TBH

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-001108  
Core Version: 8105.48  
03/05/2015 14:33:49

Standard Lot#: 25913100A2  
Last Changed By: S. CLONE #50616

DPS: A. AUTRET #53546  
TUCSON PD

Operator: A. AUTRET #53546  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:34:59
Diagnostic Test	Pass	14:35:36
Air Blank	0.000	14:36:04
0.100 Cal Check	0.097	14:36:24
Air Blank	0.000	14:36:54
Subject Test	0.257	14:37:13
Air Blank	0.000	14:37:49
Wait		
Air Blank	0.000	14:42:41
Subject Test	0.000	14:42:58
Air Blank	0.000	14:43:25
Wait		
Air Blank	0.000	14:48:26
Subject Test	XXX*	14:48:40
Air Blank	0.000	14:49:13

Run out of breath for  
Invalid sample. waited  
5 minutes, 0.000 result.

\*Invalid Sample  
Begin new deprivation period

Not a Successfully  
Completed  
Test Sequence

6 of 8  
MA  
TBH

INTOXILYZER 8100  
Location: TUCSON PD  
Serial Number: 80-001108  
Core Version: 8105.48  
03/05/2015 14:50:29

Standard Lot#: 25813100A2  
Last Changed By: S. CIONE #50616

QAS: A. AUTRET #53546  
TUCSON PD

Operator: A. AUTRET #53546  
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:50:59
0.100 Cal Check	0.098	14:51:27
Air Blank	0.000	14:51:50

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY


STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

2 of 8  
JL  
TS4  
3.11.2015

QA Specialist Amy Autret 53546 AGENCY Tucson Police Department  
DATE 3/5/15 TIME 1403  
Intoxilyzer Serial # 80-001108 Location Lab

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C + 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 097 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

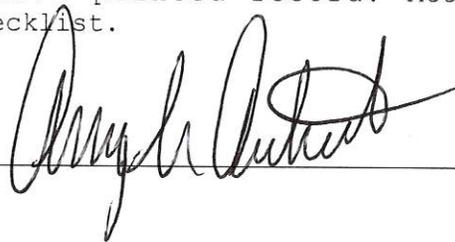SIGNATURE 

DPS Form Exh G-2 (Rev. 05-01)

8 of 8  
dd  
TS4

INTOXILYZER 8000  
Instrument Initialization  
13:37:49 03/05/2015

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-001108  
Firmware Version: 8105.48  
03/05/2015 14:03:30

Standard Lot#: 25913100A2  
Last Changed By: S. CIONE #50616

CAS: A. AUTRET #53546  
TUCSON PD

Operator: A. AUTRET #53546  
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:04:01
0.00 Cal Check	0.097	14:04:21
Air Blank	0.000	14:04:51

Preventative Maintenance  
Performed -- Timer Reset