



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001108 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

*TS4
8.6.2015*

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 7/23/2015 Cal 1 0.096 AC Cal 2 0.096 AC

Between the dates of 6/25/2015 and 7/23/2015 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on _____.
The instrument was returned to service on _____.

William L Loomis

William L. Loomis

7/27/2015

ADPS QAS Permit # 33720

SUBSCRIBED AND SWORN to before me this 27th, day of July, 2015.

Elma Escarcega

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

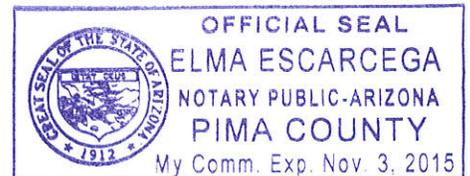


EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7
WLL
TSL 8-16-2015

QA SPECIALIST W. LOOMIS 53545 AGENCY Tucson Police Department

DATE 7/23/2015 TIME 13:37

INTOXILYZER SERIAL # 80-001108 LOCATION CRIME LAB

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 096 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE

Will Loomis

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-001108
Core Version: 8105.48
07/23/2015

13:39:59

Standard Lot#: 25913100A2
Last Changed By: S. CIONE #50616

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SQAP, TEST

DOB: 11/11/1987

Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:40:49
Diagnostic Test	Pass	13:41:24
Air Blank	0.000	13:41:52
0.100 Cal Check	0.097	13:42:13
Air Blank	0.000	13:42:42
Subject Test	0.000	13:42:57
Air Blank	0.000	13:43:25
Wait	ABT*	
Air Blank	0.000	13:43:59

*Sequence Aborted

2 of 7
WLV
104

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-001108
Core Version: 8105.48
07/23/2015

13:44:05

Standard Lot#: 25913100A2
Last Changed By: S. CIONE #50616

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SQAP,TEST

DOB: 11/11/1987

Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:45:00
Diagnostic Test	Pass	13:45:35
Air Blank	0.000	13:46:03
0.100 Cal Check	RFI*	13:46:17
Air Blank	IPS**	13:46:26

*RFI Detect

**Improper Sample

3 of 7
WLL
TSh

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-001108
Core Version: 8105.48
07/23/2015

13:46:52

Standard Lot#: 25913100A2
Last Changed By: S. CIONE #50616

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SQAP,TEST
DOB: 11/11/1987
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:47:45
Diagnostic Test	Pass	13:48:20
Air Blank	0.000	13:48:48
0.100 Cal Check	0.098	13:49:08
Air Blank	0.000	13:49:38
Subject Test	0.187*	13:52:45
Air Blank	0.000	13:53:18
Air Blank	0.000	13:53:45
Subject Test	XXX**	13:54:07
Air Blank	0.000	13:54:41

*Deficient Sample

**Invalid Sample

Begin new deprivation period

4 of 7
WV
TS4

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-001108
Core Version: 8105.48
07/23/2015 13:55:11

Standard Lot#: 25913100A2
Last Changed By: S. CIONE #50616

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	13:55:42
0.100 Cal Check	0.096	13:56:02
Air Blank	0.000	13:56:32

Preventative Maintenance
Performed -- Timer Reset

5 of 7
WLL
TS4

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7
WLL
TS4
8.6.2015

QA Specialist W. LOOMIS 53545 AGENCY Tucson Police Department
DATE 7/23/2015 TIME 13:37
Intoxilyzer Serial # 80-001108 Location CRIME LAB

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 096 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Will Loomis

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-001108
Core Version: 8105.48
07/23/2015 13:37:36

Standard Lot#: 25913100A2
Last Changed By: S. CIONE #50616

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	13:38:05
0.100 Cal Check	0.096	13:38:26
Air Blank	0.000	13:38:56

Preventative Maintenance
Performed -- Timer Reset

7 of 7
TS4
wr