



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001108 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

*TS4
9.22.2016*

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 9/15/2016 Cal 1 0.097 AC Cal 2 0.097 AC

Between the dates of 8/18/2016 and 9/15/2016 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on _____.
The instrument was returned to service on _____.

William L Loomis

William L. Loomis

9/19/2016


ADPS QAS Permit # 33720

SUBSCRIBED AND SWORN to before me this 19th day of September, 2016.

Elma Escarcega

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

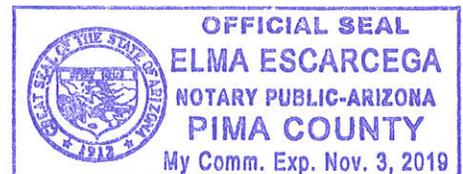

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7
WLL
BL 9.22.2016

QA SPECIALIST W. LOOMIS 53545 AGENCY Tucson Police Department
DATE 9/15/2016 TIME 08:50
INTOXILYZER SERIAL # 80-001108 LOCATION ODM

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 097 AC

Instrument is operating properly and accurately. YES NO

COMMENTS _____

SIGNATURE William Loomis

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-001108
Core Version: 8105.48
09/15/2016

08:52:33

Standard Lot#: 12415100A4
Last Changed By: W. LOOMIS #53545

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:53:30
Diagnostic Test	Pass	08:54:05
Air Blank	0.000	08:54:32
0.100 Cal Check	0.098	08:54:53
Air Blank	0.000	08:55:23
Subject Test	0.000	08:55:40
Air Blank	0.000	08:56:07
Wait	ABT*	
Air Blank	0.000	08:56:41

*Sequence Aborted

2 of 7
WU
TS4

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-001108
Core Version: 8105.48
09/15/2016

08:56:44

Standard Lot#: 12415100A4
Last Changed By: W. LOOMIS #53545

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:57:42
Diagnostic Test	Pass	08:58:17
Air Blank	0.000	08:58:45
0.100 Cal Check	RFI*	08:58:58
Air Blank	IPS**	08:59:06

*RFI Detect

**Improve Sample

3 of 7
0 WLL
TS4

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-001108
Core Version: 8105.48
09/15/2016 08:59:22

4 of 7
6 WLL
TS4

Standard Lot#: 12415100A4
Last Changed By: W. LOOMIS #53545

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:00:14
Diagnostic Test	Pass	09:00:49
Air Blank	0.000	09:01:17
0.100 Cal Check	0.098	09:01:37
Air Blank	0.000	09:02:07
Subject Test	0.084*	09:05:14
Air Blank	0.000	09:05:46
Air Blank	0.000	09:06:13
Subject Test	XXX**	09:06:36
Air Blank	0.000	09:07:09

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-001108
Core Version: 8105.48
09/15/2016 09:07:27

Standard Lot#: 12415100A4
Last Changed By: W. LOOMIS #53545

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:07:59
0.100 Cal Check	0.097	09:08:20
Air Blank	0.000	09:08:50

Preventative Maintenance
Performed -- Timer Reset

5 of 7
WLL
TS4

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7
WLU
B4 9.22.2016

QA Specialist W. LOOMIS 53545 AGENCY Tucson Police Department
DATE 9/15/2016 TIME 08:50
Intoxilyzer Serial # 80-001108 Location ODM

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C + 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 097 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Will Loomis

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-001108
Core Version: 8105.48
09/15/2016 08:50:24

Standard Lot#: 12415100A4
Last Changed By: W. LOOMIS #53545

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	08:50:54
0.100 Cal Check	0.097	08:51:15
Air Blank	0.000	08:51:45

Preventative Maintenance
Performed -- Timer Reset

7 of 7
0 WLU
104