



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

ju

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001163 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 12-1-2006 Cal 1 0.099 AC Cal 2 0.099 AC


The instrument was placed into service on 12-1-2006.

Terry S. Gallegos
Terry S. Gallegos

1-4-2007

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 4, day of January, 2007.

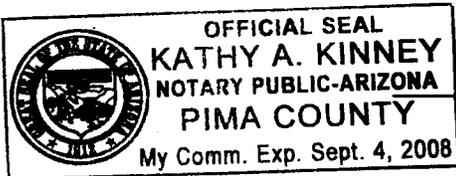

Kathy A. Kinney
Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

1029
DJ

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

jk

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 12-1-2006 TIME 1559

INTOXILYZER SERIAL # 80-001163 LOCATION Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0.000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed ✓✓
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed ✓
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0.100 AC
Result: 0.099 AC

Instrument is operating properly and accurately. YES NO

COMMENTS During deficient sample check, I introduced the
breath sample too long and it was accepted as a valid test.

SIGNATURE 

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

20 of 9
B4

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001163
12/01/2006

16:00:31

Standard Lot#: 502602I
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:01:23
Diagnostic Test	Pass	16:01:51
Air Blank	0.000	16:02:18
0.100 Cal Check	0.099	16:02:37
Air Blank	0.000	16:03:06
Subject Test	0.000	16:03:29
Air Blank	0.000	16:03:56
Five Minute Wait	ABT*	
Air Blank	0.000	16:05:20

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

3029
TJ

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001163
12/01/2006 16:06:26

Standard Lot#: 502602I
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP,TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	IPS*	16:07:25
Air Blank	0.000	16:07:53

*Improper Sample

Not a Successfully Completed Test Sequence

4 of 9
TS4

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001163
12/01/2006

16:07:59

Standard Lot#: 502602I
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP,TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:08:49
Diagnostic Test	Pass	16:09:16
Air Blank	0.000	16:09:43
0.100 Cal Check	0.099	16:10:02
Air Blank	0.000	16:10:32
Subject Test	RFI*	16:10:43
Air Blank	0.000	16:11:11

*RFI Detect

Not a Successfully Completed Test Sequence

5 of 9
Tsh

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001163
12/01/2006

16:11:18

Standard Lot#: 502602I
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:12:12
Diagnostic Test	Pass	16:12:39
Air Blank	0.000	16:13:07
0.100 Cal Check	0.099	16:13:26
Air Blank	0.000	16:13:55
Subject Test	0.017	16:14:16
Air Blank	0.000	16:14:47
Five Minute Wait	ABT*	
Air Blank	0.000	16:15:33

*Sequence Aborted

Not a Successfully Completed Test Sequence

6 of 9
PH

INTOXILYZER 8000 *je*
Location: ARIZONA
Serial Number : 80-001163
12/01/2006

16:16:02

Standard Lot#: 502602I
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP,TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:16:53
Diagnostic Test	Pass	16:17:20
Air Blank	0.000	16:17:47
0.100 Cal Check	0.099	16:18:06
Air Blank	0.000	16:18:36
Subject Test	0.174*	16:21:44
Air Blank	0.000	16:22:16
Air Blank	0.000	16:22:43
Subject Test	XXX**	16:24:04
Air Blank	0.000	16:24:39

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

70f9
TSh

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001163
12/01/2006

16:25:49

Standard Lot#: 502602I
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	16:26:16
Cal Check	0.099	16:26:35
Air Blank	0.000	16:27:05


Operator's Signature

8 of 9
TDH

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD CALIBRATION CHECK PROCEDURE

je

QA Specialist Terry Gallegos 28047 AGENCY Tucson Police Department
DATE 12-1-2006 TIME 1559
Intoxilyzer Serial # 80-001163 Location Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE 

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

9 of 9
Pg

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001183
12/01/2006 15:58:34

Standard Lot#: 502602I
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	15:59:02
Cal Check	0.099	15:59:21
Air Blank	0.000	15:59:50


Operator's Signature



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

INTOXILYZER SERVICE RECORD

ju
1 of 1
TJ4

Intoxilyzer Serial Number: 80-001163

Initial Service Date: 12.1.2006

CMI, Inc. Certificate of Calibration Date: 2.4.2005