



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory

WMM
6-5-07

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001163 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 5-24-2007 Cal 1 0.100 AC Cal 2 0.098 AC

Between the dates of 4-26-2007 and 5-24-2007 the instrument was operating properly and accurately.


The instrument was taken out of service for repair on .
The instrument was returned to service on .

Terry S. Gallegos
Terry S. Gallegos

6.5.2007

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 5, day of June, 2007.

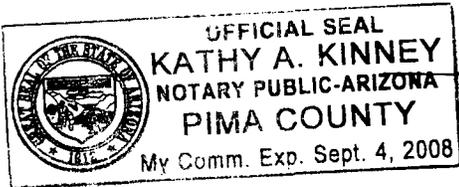

Kathy A Kinney
Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 8
104 mm
6-5-03

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 5-24-2007 TIME 1735

INTOXILYZER SERIAL # 80-001163 LOCATION DUI HQ

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed ✓✓
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 098 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Terry Gallegos

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001163
05/24/2007 17:40:46

2 of 8
TS4
6-5-07

Standard Lot#: 635301I
Last Changed By: G. EPPLEY #26227

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|------------------|--------|----------|
| Air Blank | 0.000 | 17:41:41 |
| Diagnostic Test | Pass | 17:42:09 |
| Air Blank | 0.000 | 17:42:36 |
| 0.100 Cal Check | 0.099 | 17:42:55 |
| Air Blank | 0.000 | 17:43:24 |
| Subject Test | 0.000 | 17:43:54 |
| Air Blank | 0.000 | 17:44:21 |
| Five Minute Wait | ABT* | |
| Air Blank | 0.000 | 17:45:24 |

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001163
05/24/2007 17:45:45

3 of 8
T4
TUCSON

Standard Lot#: 635301I
Last Changed By: G. EPPLEY #26227

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 17:46:40 |
| Diagnostic Test | Pass | 17:47:07 |
| Air Blank | 0.000 | 17:47:35 |
| 0.100 Cal Check | 0.100 | 17:47:53 |
| Air Blank | RFI* | 17:48:00 |
| Air Blank | IPS** | 17:48:09 |

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001163
05/24/2007 17:48:30

4 of 8
TSy 2ms
6-5-07

Standard Lot#: 635301I
Last Changed By: G. EPPLEY #26227

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------|--------|----------|
| Air Blank | AMB* | 17:49:18 |
| Air Blank | 0.000 | 17:49:46 |

*Ambient Fail

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-001163
05/24/2007 17:49:56

5 of 8
Dy
6-5-07

Standard Lot#: 635301I
Last Changed By: G. EPPLEY #26227

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 17:50:47 |
| Diagnostic Test | Pass | 17:51:14 |
| Air Blank | 0.000 | 17:51:41 |
| 0.100 Cal Check | 0.099 | 17:52:00 |
| Air Blank | 0.000 | 17:52:30 |
| Subject Test | 0.129* | 17:55:38 |
| Air Blank | 0.000 | 17:56:11 |
| Air Blank | 0.000 | 17:56:38 |
| Subject Test | XXX** | 17:57:02 |
| Air Blank | 0.000 | 17:57:36 |

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001163
05/24/2007 17:58:00

6 of 8
TSY
with
500

Standard Lot#: 635301I
Last Changed By: G. EPPLEY #26227

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

| Test | g/210L | Time |
|-----------|--------|----------|
| Air Blank | 0.000 | 17:58:27 |
| Cal Check | 0.098 | 17:58:46 |
| Air Blank | 0.000 | 17:59:15 |



Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

*7 of 8
TSG
6-3-07*

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 5.24.2007 TIME 1735
Intoxilyzer Serial # 80-001163 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE *Terry Gallegos*

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department


INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001163
05/24/2007 17:35:46

8 of 8
TS4
6/5/07

Standard Lot#: 635301I
Last Changed By: G. EPPLEY #26227

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

| Test | g/210L | Time |
|-----------|--------|----------|
| Air Blank | 0.000 | 17:36:14 |
| Cal Check | 0.100 | 17:36:32 |
| Air Blank | 0.000 | 17:37:02 |

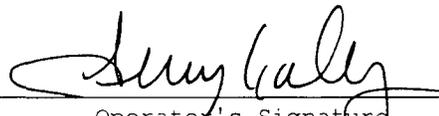

Operator's Signature