



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001163 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4
8-8-07

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 7-31-07 Cal 1 0.098 AC Cal 2 0.098 AC

Between the dates of 6-21-07 and 7-31-07 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.


Clare Conley

Clare L. Conley

8-3-07

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 3, day of August, 2007.

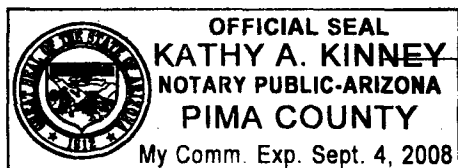

Kathy A Kinney
Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

TS4
8-807

QA SPECIALIST Clare Conley AGENCY Tucson Police Department
DATE 7-31-07 TIME 0824
INTOXILYZER SERIAL # 80-001163 LOCATION TPD Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
(✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 098 AC

Instrument is operating properly and accurately. YES ☒ NO ☐

COMMENTS Instrument not available prior to this date.

SIGNATURE Clare Conley

2 of 8 cu

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-001163

07/31/2007

08:26:35

ph

Standard Lot#: 635301I

Last Changed By: G. EPPLEY

#26227

QAS: C. CONLEY

#34517

TUCSON PD

Operator: C. CONLEY

#34517

TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:27:36
Diagnostic Test	Pass	08:28:03
Air Blank	0.000	08:28:31
0.100 Cal Check	0.099	08:28:49
Air Blank	0.000	08:29:19
Subject Test	0.000	08:29:34
Air Blank	ABT*	08:29:44
Air Blank	0.000	08:30:13

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

348u

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-001163

07/31/2007

08:33:16

PM

Standard Lot#: 635301I

Last Changed By: G. EPPLEY

#26227

QAS: C. CONLEY

#34517

TUCSON PD

Operator: C. CONLEY

#34517

TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:34:14
Diagnostic Test	Pass	08:34:41
Air Blank	0.000	08:35:09
0.100 Cal Check	RFI*	08:35:20
Air Blank	0.000	08:35:51

*RFI Detect

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

4 of 8 Cu

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-001163

07/31/2007

08:36:11

15h

Standard Lot#: 635301I

Last Changed By: G. EPPLEY

#26227

QAS: C. CONLEY

#34517

TUCSON PD

Operator: C. CONLEY

#34517

TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	IPS*	08:36:46
Air Blank	0.000	08:37:15

*Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

5 of 8 cu

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-001163

07/31/2007

08:38:05

10h

Standard Lot#: 635301I

Last Changed By: G. EPPLEY

#26227

QAS: C. CONLEY

#34517

TUCSON PD

Operator: C. CONLEY

#34517

TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:38:59
Diagnostic Test	Pass	08:39:26
Air Blank	0.000	08:39:54
0.100 Cal Check	0.098	08:40:13
Air Blank	0.000	08:40:42
Subject Test	0.063*	08:43:51
Air Blank	0.000	08:44:23
Air Blank	0.000	08:44:50
Subject Test	XXX**	08:45:08
Air Blank	0.000	08:45:42

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

6 of 8 cu

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001163

07/31/2007

08:46:12

104

Standard Lot#: 635301I

Last Changed By: G. EPPLEY #26227

QAS: C. CONLEY #34517

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	08:46:39
Cal Check	0.098	08:46:58
Air Blank	0.000	08:47:27

Clare Conley

Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

7 of 8 cu

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

B4

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Clare Conley AGENCY Tucson Police Department
DATE 7-31-07 TIME 0824
Intoxilyzer Serial # 80-001163 Location TPD Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$.
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 098 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

8 of 8 cu

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001163

07/31/2007

08:24:21

ph

Standard Lot#: 635301I

Last Changed By: G. EPPLEY #26227

QAS: C. CONLEY #34517

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	08:24:48
Cal Check	0.098	08:25:07
Air Blank	0.000	08:25:36

Clare Conley
Operator's Signature