



CERTIFICATION by CUSTODIAN of RECORDS  
Tucson Police Department Crime Laboratory

8.31.07

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-001163 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 8-17-2007 Cal 1 0.097 AC Cal 2 0.097 AC

Between the dates of 7.31.2007 and 8.17.2007 the instrument was operating properly and accurately.


The instrument was taken out of service for repair on \_\_\_\_\_.  
The instrument was returned to service on \_\_\_\_\_.

Terry S. Gallegos  
Terry S. Gallegos

8.30.2007

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 30, day of August, 2007.

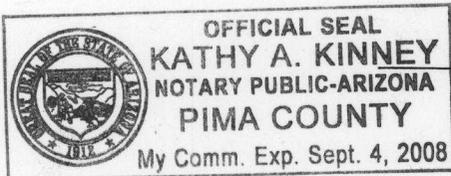

Kathy A Kinney  
Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

CU 8-31-07

1 of 7  
134

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 8.17.2007 TIME 1555

INTOXILYZER SERIAL # 80-001163 LOCATION DNI HQ

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed ✓✓
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 097 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---




---



---



---

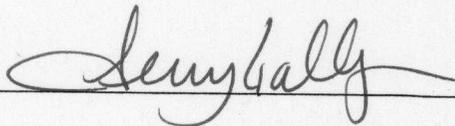SIGNATURE 

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

Oh

20f7  
B4

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-001163  
08/17/2007

15:59:11

Standard Lot#: 706601I  
Last Changed By: S. CLEARY #44524

QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:59:59
Diagnostic Test	Pass	16:00:27
Air Blank	0.000	16:00:54
0.100 Cal Check	0.097	16:01:13
Air Blank	0.000	16:01:42
Subject Test	0.000	16:02:02
Air Blank	0.000	16:02:29
Five Minute Wait	ABT*	
Air Blank	0.000	16:03:16

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-001163  
08/17/2007

16:03:22

*CM*

*30f7  
TS4*

---

Standard Lot#: 706601I  
Last Changed By: S. CLEARY #44524

---

QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

---

Operator: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

---

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	16:04:16
Diagnostic Test	Pass	16:04:45
Air Blank	0.000	16:05:12
0.100 Cal Check	0.097	16:05:31
Air Blank	0.000	16:06:01
Subject Test	RFI*	16:06:12
Air Blank	IPS**	16:06:26

\*RFI Detect

\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-001163  
08/17/2007

Oh

16:09:38

4 of 7  
TS4

---

Standard Lot#: 706601I  
Last Changed By: S. CLEARY #44524

---

QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

---

Operator: T. GALLEGOS #28847  
TUCSON PD CRIME LAB

---

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	16:10:28
Diagnostic Test	Pass	16:10:55
Air Blank	0.000	16:11:22
0.100 Cal Check	0.098	16:11:41
Air Blank	0.000	16:12:11
Subject Test	0.344*	16:15:19
Air Blank	0.000	16:15:52
Air Blank	0.000	16:16:19
Subject Test	XXX**	16:17:07
Air Blank	0.000	16:17:41

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001163  
08/17/2007

Oh

16:18:49

-----  
Standard Lot#: 706601I  
Last Changed By: S. CLEARY #44524  
-----

507  
734

-----  
QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB  
-----

Test	g/210L	Time
Air Blank	0.000	16:19:16
Cal Check	0.097	16:19:35
Air Blank	0.000	16:20:05

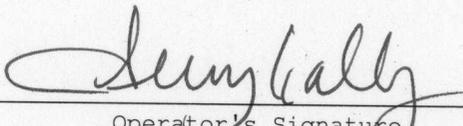
  
-----  
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A) *On 8-31-07*

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

*6 of 7  
TSK*

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 28897 AGENCY Tucson Police Department

DATE 8-17-2007 TIME 1555

Intoxilyzer Serial # 80-001163 Location DUI HQ

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 097 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE *Terry Gallegos*

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001163  
08/17/2007 15:55:54

CM

7077  
74

-----  
Standard Lot#: 706601I  
Last Changed By: S. CLEARY #44524  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD CRIME LAB  
-----

Test	g/210L	Time
Air Blank	0.000	15:56:21
Cal Check	0.097	15:56:40
Air Blank	0.000	15:57:09

  
-----  
Operator's Signature