



CERTIFICATION by CUSTODIAN of RECORDS  
Tucson Police Department Crime Laboratory

CM  
2-8-08

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-001163 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 2/8/08 Cal 1 0.099 AC Cal 2 0.099 AC

Between the dates of 1/10/08 and 2/8/08 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
The instrument was returned to service on -.


Sarah Cione

Sarah L. Cione

2/8/08

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 8, day of February, 2008.

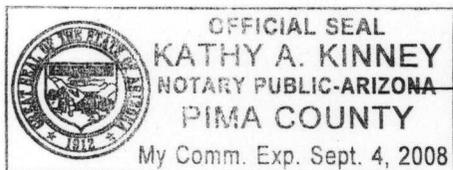

Kathy A Kinney

Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

1 of 75c

ARIZONA DEPARTMENT OF PUBLIC SAFETY

DM

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

2-8-08

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 2/8/08 TIME 10:11

INTOXILYZER SERIAL # 80-00163 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 099 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE

Sarah Cione

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-001163  
02/08/2008

10:13:01

2 of 7sc  
CU  
2-8-08

---

Standard Lot#: 719802I  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP,TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	10:13:50
Diagnostic Test	Pass	10:14:17
Air Blank	0.000	10:14:45
0.100 Cal Check	0.100	10:15:04
Air Blank	0.000	10:15:33
Subject Test	0.000	10:16:02
Air Blank	0.000	10:16:29
Five Minute Wait	ABT*	
Air Blank	0.000	10:17:15

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-001163  
02/08/2008

10:17:28

3 of 75C

OH

2-8-08

---

Standard Lot#: 719802I  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	RFI*	10:18:06
Air Blank	IPS**	10:18:14

\*RFI Detect

\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-001163  
02/08/2008

10:19:03

4 of 75c

CU

2-8-08

---

Standard Lot#: 719802I  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, TEST  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	10:19:52
Diagnostic Test	Pass	10:20:19
Air Blank	0.000	10:20:47
0.100 Cal Check	0.100	10:21:05
Air Blank	0.000	10:21:35
Subject Test	0.061*	10:24:43
Air Blank	0.000	10:25:16
Air Blank	0.000	10:25:43
Subject Test	XXX**	10:26:06
Air Blank	0.000	10:26:40

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001163  
02/08/2008 10:29:00

5 of 7sc

lu

2-8-08

-----  
Standard Lot#: 719802I  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	10:29:27
Cal Check	0.099	10:29:46
Air Blank	0.000	10:30:15

*Sarah Cione*

Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

6 of 75c

CM  
2-8-08

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

## STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 2/8/08 TIME 10:11  
Intoxilyzer Serial # 80-001163 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach  $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001163  
02/08/2008 10:11:37

7 of 7sc

du  
2-8-08

-----  
Standard Lot#: 719802I  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	10:12:05
Cal Check	0.099	10:12:23
Air Blank	0.000	10:12:53

*Sarah Cione*

\_\_\_\_\_  
Operator's Signature



TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

1 of 2 SC

lu  
2-8-08

INTOXILYZER SERVICE RECORD  
FOR  
CALIBRATION GAS TANK CHANGE

Date: 2/6/08  
 QAS: Sarah Cione  
 Intoxilyzer Serial Number: 80-001163  
 New Lot Number: 719802 I  
 New Expiration Date: 7/23/2009

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."  
 Select "D" for dry gas and enter the following information:  
     Enter Target Value: 0.100  
     Enter Control Lot #: From the tank label  
     Enter Exp Date: From the tank label (mm/dd/yyyy)  
     Edit info Y/N? Y or N  
     Update Std Info? Y ENTER  
     "Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing."  
 Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 100 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove and sign printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

2 of 2 SC


02-8-08

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001163  
02/06/2008 16:53:28

-----  
Standard Lot#: 719802I  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	16:53:55
Cal Check	0.100	16:54:14
Air Blank	0.000	16:54:43

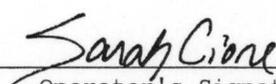  
-----  
Operator's Signature