



CERTIFICATION by CUSTODIAN of RECORDS

*Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917*

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001163 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

184
6-10-09

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 06/04/09 Cal 1 0.100 AC Cal 2 0.100 AC

Between the dates of 05/07/09 and 06/04/09 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.
The instrument was returned to service on —.

Cheryl K. Straub

Cheryl K. Straub


06/09/09

ADPS QAS Permit # 33408

SUBSCRIBED AND SWORN to before me this 9th, day of June, 2009.

A.M. Hernández

Notary Public

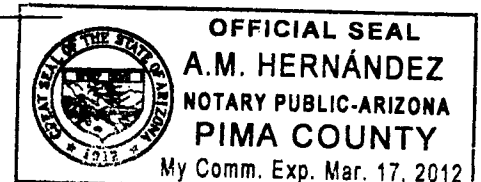

*Includes documents specified in A.A.C. R13-10-104.A.5

177
CS

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

TSY
6-10-09

QA SPECIALIST Cheryl STRAUB AGENCY Tucson Police Department
DATE 06/04/09 TIME 08:47
INTOXILYZER SERIAL # 80-001163 LOCATION DUI #2

(☒) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (☒) 1. Clock time check.
(☒) 2. Date check.

OPERATIONAL TESTS

- (☒) 1. Alcohol-free subject test result 0. 000 AC.
(☒) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
(☒) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(☒) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 100 AC

Instrument is operating properly and accurately. YES X NO

COMMENTS

SIGNATURE

Cheryl Straub

277
c8

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001163
06/04/2009

09:04:46

134

Standard Lot#: 804302I
Last Changed By: C. CONLEY #34517

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: TEST,TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:05:41
Diagnostic Test	Pass	09:06:09
Air Blank	0.000	09:06:36
0.100 Cal Check	0.102	09:06:55
Air Blank	0.000	09:07:25
Subject Test	0.000	09:07:43
Air Blank	ABT*	09:07:50
Air Blank	0.000	09:08:18

*Sequence Aborted

Not a Successfully Completed Test Sequence

387
C8

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001163
06/04/2009

08:52:34

154

Standard Lot#: 804302I
Last Changed By: C. CONLEY #34517

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: TEST,TEST

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	08:53:20
Air Blank	IPS**	08:53:27

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1527, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

487
CS

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001163
06/04/2009

09:09:28

AS4

Standard Lot#: 804302I
Last Changed By: C. CONLEY #34517

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: TEST,TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:10:32
Diagnostic Test	Pass	09:10:59
Air Blank	0.000	09:11:27
0.100 Cal Check	0.101	09:11:45
Air Blank	0.000	09:12:15
Subject Test	0.000*	09:15:24
Air Blank	0.000	09:15:51
Air Blank	0.000	09:16:18
Subject Test	XXX**	09:16:35
Air Blank	0.000	09:17:08

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1527, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department

587
CS

INTOXILYZER 8000

Location: ARIZONA

Serial Number: 80-001163

06/04/2009

09:19:12

154

Standard Lot#: 804302I

Last Changed By: C. CONLEY #34517

QAS: C. STRAUB #50772

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:19:40
Cal Check	0.100	09:19:58
Air Blank	0.000	09:20:28



Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Cheryl Straub AGENCY Tucson Police Department
DATE 06/04/09 TIME 08:47
Intoxilyzer Serial # 80-001163 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution 100 AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Cheryl Straub

DPS Form Exh G-2 (Rev. 05-01)

787
CS


INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001163
06/04/2009 08:47:50

154

Standard Lot#: 804302I
Last Changed By: C. CONLEY #34517

QAS: C. STRAUB #50772
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	08:48:18
Cal Check	0.100	08:48:36
Air Blank	0.000	08:49:06

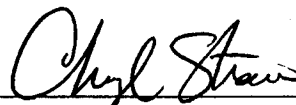

Operator's Signature