



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001163 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4
6.26.09

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6/25/09 Cal 1 0.100 AC Cal 2 0.099 AC

Between the dates of 6/4/09 and 6/25/09 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -
The instrument was returned to service on -

Sarah Cione

Sarah L. Cione

6/25/09


ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 25th, day of June, 20 09.

A.M. Hernandez

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

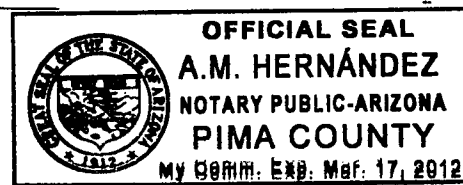

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 7 SC

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

By 6/26/09

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department
DATE 6/25/09 TIME 8:38
INTOXILYZER SERIAL # 80-001163 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
(✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
() 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES ✓ NO

COMMENTS

Tank was not changed, just accidentally updated standard
information (pg. 5 of 7)

SIGNATURE

Sarah Cione

2 of 7sc

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001163
06/25/2009

08:40:23

B4

Standard Lot#: 804302I
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP,TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:41:18
Diagnostic Test	Pass	08:41:46
Air Blank	0.000	08:42:13
0.100 Cal Check	0.100	08:42:32
Air Blank	0.000	08:43:02
Subject Test	0.000	08:43:41
Air Blank	0.000	08:44:09
Five Minute Wait	ABT*	
Air Blank	0.000	08:44:54

*Sequence Aborted

Not a Successfully Completed Test Sequence

3 of 7_{sc}

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001163
06/25/2009

08:45:01

Standard Lot#: 804302I
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP,TEST

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	08:45:59
Air Blank	IPS**	08:46:16

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

4 of 7_{sc}

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001163
06/25/2009

08:46:49

by

Standard Lot#: 804302I
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP,TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:47:51
Diagnostic Test	Pass	08:48:19
Air Blank	0.000	08:48:46
0.100 Cal Check	0.100	08:49:05
Air Blank	0.000	08:49:35
Subject Test	0.074*	08:52:43
Air Blank	0.000	08:53:15
Air Blank	0.000	08:53:43
Subject Test	XXX**	08:54:04
Air Blank	0.000	08:54:37

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

5 of 7sc

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001163
06/25/2009 08:55:33

Standard Lot#: 804302I
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	08:56:00
Cal Check	0.099	08:56:19
Air Blank	0.000	08:56:49


Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

6 of 7_{sc}

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000TSH
6-26-09

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 6/25/09 TIME 8:38
Intoxilyzer Serial # 80-001163 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: ARIZONA

Serial Number: 80-001163

06/25/2009

08:38:50

Standard Lot#: 804302ILast Changed By: C. CONLEY #34517

QAS: S. CIONE #50616

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	08:39:18
Cal Check	0.100	08:39:36
Air Blank	0.000	08:40:06



Operator's Signature