



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001143 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 07/23/09 Cal 1 0.100 AC Cal 2 0.101 AC

Between the dates of 06/25/09 and 07/23/09 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on .
The instrument was returned to service on .

Cheryl K. Straub

Cheryl K. Straub

07/28/09


ADPS QAS Permit # 33408

SUBSCRIBED AND SWORN to before me this 28, day of July, 2009.

A.M. Hernandez

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

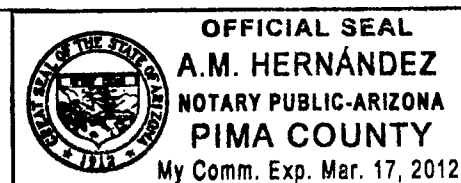

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

154
7-28-09

QA SPECIALIST Cheryl STRAUB AGENCY Tucson Police Department

DATE 07/23/09 TIME 09:08

INTOXILYZER SERIAL # 80-001163 LOCATION TEAM II/DUI

(☒) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (☒) 1. Clock time check.
(☒) 2. Date check.

OPERATIONAL TESTS

- (☒) 1. Alcohol-free subject test result 0. 000 AC.
(☒) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
(☒) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(☒) 4. Standard Calibration Check standard 0. 160 AC
Result: 0. 101 AC

Instrument is operating properly and accurately. YES ☒ NO ☐

COMMENTS

SIGNATURE

Cheryl Straub

DPS Form Exh G-4 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

218cs

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001163
07/23/2009

09:10:52

134

Standard Lot#: 804302I
Last Changed By: S. CIONE #50616

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:11:52
Diagnostic Test	Pass	09:12:19
Air Blank	0.000	09:12:46
0.100 Cal Check	0.100	09:13:05
Air Blank	0.000	09:13:35
Subject Test	0.000	09:14:29
Air Blank	ABT*	09:14:36
Air Blank	0.000	09:15:05

*Sequence Aborted

Not a Successfully Completed Test Sequence

398cs

154

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001163
07/23/2009

09:15:43

Standard Lot#: 804302I
Last Changed By: S. CIONE #50616

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	09:16:30
Air Blank	IPS**	09:16:49

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

498 CS

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001163
07/23/2009

09:17:27

104

Standard Lot#: 804302I
Last Changed By: S. CIONE #50616

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:18:24
Diagnostic Test	Pass	09:18:52
Air Blank	0.000	09:19:19
0.100 Cal Check	0.100	09:19:38
Air Blank	0.000	09:20:08
Subject Test	0.039*	09:23:16
Air Blank	0.000	09:23:47
Air Blank	0.000	09:24:15
Subject Test	0.181	09:24:56
Air Blank	0.000	09:25:30
Five Minute Wait	ABT**	
Air Blank	0.000	09:27:28

*Deficient Sample

**Sequence Aborted

*while performing Deficient sample test I
did not have enough breath spray in my
mouth and blew too long, instrument took it
As a Subject Test. Test repeated.
CS*

Not a Successfully Completed Test Sequence

57805

54

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001163
07/23/2009

09:27:47

Standard Lot#: 804302I
Last Changed By: S. CIONE #50616

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:28:51
Diagnostic Test	Pass	09:29:18
Air Blank	0.000	09:29:46
0.100 Cal Check	0.100	09:30:04
Air Blank	0.000	09:30:34
Subject Test	XXX*	09:31:04
Air Blank	0.000	09:31:38

*Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: ARIZONA

Serial Number: 80-001163

07/23/2009

09:50:18


PS4

Standard Lot#: 804302ILast Changed By: S. CIONE #50616

QAS: C. STRAUB #50772

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:50:45
Cal Check	0.101	09:51:04
Air Blank	0.000	09:51:34

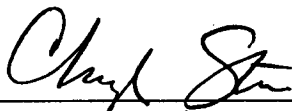_____
Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

7/28/09

QA Specialist Cheryl Straub AGENCY Tucson Police Department
DATE 07/23/09 TIME 09:08
Intoxilyzer Serial # 80-001163 Location TEAM II/DUT

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Cheryl Straub

DPS Form Exh G-2 (Rev. 05-01)

8 g 8 cs


INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001163
07/23/2009 09:08:23

B4

Standard Lot#: 804302I
Last Changed By: S. CIONE #50616

QAS: C. STRAUB #50772
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:08:50
Cal Check	0.100	09:09:09
Air Blank	0.000	09:09:39

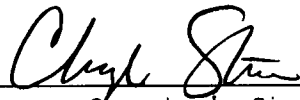

Operator's Signature