



# CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-001163 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4  
7-6-2010

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 07/01/10 Cal 1 0.099 AC Cal 2 0.099 AC

Between the dates of 06/03/10 and 07/01/10  
the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.  
The instrument was returned to service on N/A.

Cheryl K. Straub

Cheryl K. Straub

07/01/10

ADPS QAS Permit # 33408

SUBSCRIBED AND SWORN to before me this 1, day of

July, 2010.

A.M. Hernandez

Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

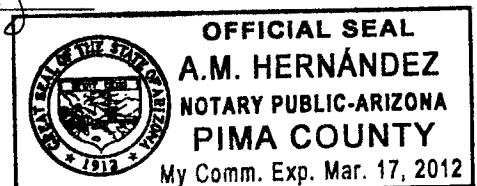


EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Cheryl STRAUB AGENCY Tucson Police Department

DATE 07/01/10 TIME 12:22

INTOXILYZER SERIAL # 80-001163 LOCATION CRIME LAB

( ☒ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ☒ ) 1. Clock time check.  
( ☒ ) 2. Date check.

OPERATIONAL TESTS

- ( ☒ ) 1. Alcohol-free subject test result 0. 000 AC.  
( ☒ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed  
( ☒ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.  
( ☒ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 099 AC

Instrument is operating properly and accurately. YES ☒ NO ☐

COMMENTS

SIGNATURE Cheryl Straub

277  
CS

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-001163  
07/01/2010

12:32:19

BM

Standard Lot#: 719919  
Last Changed By: S. CIONE #50616

QAS: C. STRAUB #50772  
TUCSON PD

Operator: C. STRAUB #50772  
TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:33:17
Diagnostic Test	Pass	12:33:44
Air Blank	0.000	12:34:12
0.100 Cal Check	0.100	12:34:30
Air Blank	0.000	12:35:00
Subject Test	0.000	12:35:26
Air Blank	ABT*	12:35:32
Air Blank	0.000	12:36:01

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

377  
CS

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-001163  
07/01/2010

12:37:12

BM

---

Standard Lot#: 719919  
Last Changed By: S. CIONE #50616

---

QAS: C. STRAUB #50772  
TUCSON PD

---

Operator: C. STRAUB #50772  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

---

Test	g/210L	Time
Air Blank	RFI*	12:37:50
Air Blank	IPS**	12:38:00

---

\*RFI Detect

\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-001163  
07/01/2010

12:38:08

Standard Lot#: 719919  
Last Changed By: S. CIONE #50616

QAS: C. STRAUB #50772  
TUCSON PD

Operator: C. STRAUB #50772  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:38:58
Diagnostic Test	Pass	12:39:26
Air Blank	0.000	12:39:53
0.100 Cal Check	0.099	12:40:12
Air Blank	0.000	12:40:42
Subject Test	0.083*	12:43:50
Air Blank	0.000	12:44:22
Air Blank	0.000	12:44:49
Subject Test	XXX**	12:45:26
Air Blank	0.000	12:45:59

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

587  
CS

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001163  
07/01/2010 12:46:15

BA

-----  
Standard Lot#: 719919  
Last Changed By: S. CIONE #50616  
-----

QAS: C. STRAUB #50772  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	12:46:42
Cal Check	0.099	12:47:01
Air Blank	0.000	12:47:31

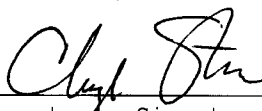
  
\_\_\_\_\_  
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Cheryl STRAUB AGENCY Tucson Police Department  
DATE 07/01/10 TIME 12:22  
Intoxilyzer Serial # 80-001163 Location CRIME LAB

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution N/A AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Cheryl Straub

DPS Form Exh G-2 (Rev. 05-01)

777  
CS

BH


1

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001163  
07/01/2010 12:22:42

Standard Lot#: 719919  
Last Changed By: S. CIONE #50616

OAS: C. STRAUB #50772  
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	12:23:09
Cal Check	0.099	12:23:28
Air Blank	0.000	12:23:58

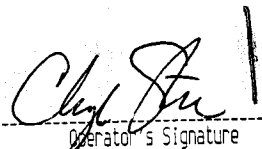  
Operator's Signature



TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

18 of 2 SC

INTOXILYZER SERVICE RECORD  
FOR  
CALIBRATION GAS TANK CHANGE

By 7-6-2010

Date: 6/23/2010

QAS: Sarah Cione

Intoxilyzer Serial Number: 80-001163

New Lot Number: 719919

New Expiration Date: 2/2/2012

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change." Select "D" for dry gas and enter the following information:  
Enter Target Value: 0.100  
Enter Control Lot #: From the tank label  
Enter Exp Date: From the tank label (mm/dd/yyyy)  
Edit info Y/N? Y or N  
Update Std Info? Y ENTER  
"Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "Control Testing." Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 100 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

INTOXILYZER 8000

Location: ARIZONA

Serial Number: 80-001163

06/23/2010

12:58:37

-----  
Standard Lot#: 719919

Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	12:59:05
Cal Check	0.100	12:59:24
Air Blank	0.000	12:59:53

Sarah Cione  
Operator's Signature