



CERTIFICATION by CUSTODIAN of RECORDS  
 Tucson Police Department Crime Laboratory  
 270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-001163 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

*134*  
*2-10-2011*

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 02/10/2011 Cal 1 0.100 AC Cal 2 0.098 AC

Between the dates of 01/13/2011 and 02/10/2011 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.  
 The instrument was returned to service on N/A.

*Cheryl K. Straub*

Cheryl K. Straub


02/10/2011

ADPS QAS Permit # 33408

SUBSCRIBED AND SWORN to before me this 10, day of February, 2011.

*A.M. Hernandez*

Notary Public

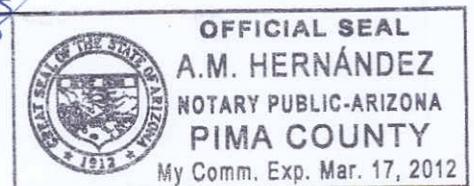

\*Includes documents specified in A.A.C. R13-10-104.A.5

167  
CS

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

TSA  
2-10-2011

QA SPECIALIST Cheryl STRAUB AGENCY Tucson Police Department

DATE 02/10/2011 TIME 06:47

INTOXILYZER SERIAL # 80-001163 LOCATION TEAM 3

(  ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (  ) 1. Clock time check.
- (  ) 2. Date check.

OPERATIONAL TESTS

- (  ) 1. Alcohol-free subject test result 0. 000 AC.
- (  ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (  ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (  ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 098 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Cheryl Straub

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

277  
CS

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-001163  
02/10/2011

06:48:55

DM

Standard Lot#: 719919  
Last Changed By: S. CIONE #50616

QAS: C. STRAUB #50772  
TUCSON PD

Operator: C. STRAUB #50772  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	06:49:44
Diagnostic Test	Pass	06:50:12
Air Blank	0.000	06:50:39
0.100 Cal Check	0.099	06:50:58
Air Blank	0.000	06:51:28
Subject Test	0.000	06:51:50
Air Blank	ABT*	06:51:56
Air Blank	0.000	06:52:25

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

327  
CS

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-001163  
02/10/2011

06:52:28

TS4

Standard Lot#: 719919  
Last Changed By: S. CIONE #50616

QAS: C. STRAUB #50772  
TUCSON PD

Operator: C. STRAUB #50772  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	06:53:00
Air Blank	IPS**	06:53:09

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

477  
CS

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-001163  
02/10/2011

06:53:17

PH

Standard Lot#: 719919  
Last Changed By: S. CIONE #50616

OAS: C. STRAUB #50772  
TUCSON PD

Operator: C. STRAUB #50772  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	06:54:05
Diagnostic Test	Pass	06:54:33
Air Blank	0.000	06:55:00
0.100 Cal Check	0.099	06:55:19
Air Blank	0.000	06:55:49
Subject Test	0.160*	06:58:57
Air Blank	0.000	06:59:29
Air Blank	0.000	06:59:57
Subject Test	XXX**	07:00:15
Air Blank	0.000	07:00:50

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

507  
CS


INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001163  
02/10/2011 07:01:08

-----  
Standard Lot#: 719919  
Last Changed By: S. CIONE #50616  
-----

TS4

-----  
QAS: C. STRAUB #50772  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	07:01:35
Cal Check	0.098	07:01:54
Air Blank	0.000	07:02:24

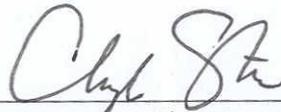

Operator's Signature

687  
CS

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY


STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

TSY  
2-10-2011

QA Specialist Cheryl STRAUB AGENCY Tucson Police Department  
DATE 02/10/2011 TIME 06:47  
Intoxilyzer Serial # 80-001163 Location TEAM 3

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution N/A AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 100 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

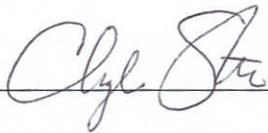SIGNATURE 

DPS Form Exh G-2 (Rev. 05-01)

777  
CS


INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001163  
02/10/2011 06:47:06

104

-----  
Standard Lot#: 719919  
Last Changed By: S. CIONE #50616  
-----

QAS: C. STRAUB #50772  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	06:47:33
Cal Check	0.100	06:47:52
Air Blank	0.000	06:48:22

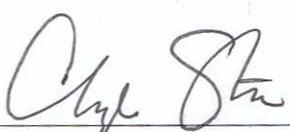

Operator's Signature