



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001163 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4
5.10.2011

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 05/05/11 Cal 1 0.099 AC Cal 2 0.098 AC

Between the dates of 04/07/2011 and 05/05/2011 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.


Cheryl K. Straub

Cheryl K. Straub

05/10/2011

ADPS QAS Permit # 33408

SUBSCRIBED AND SWORN to before me this 10, day of May, 20 11.

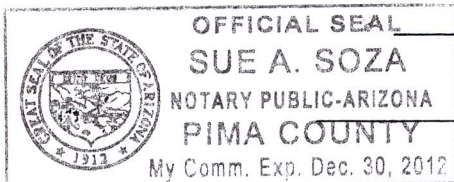

Sue A. Soza

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

187
CS

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

BS4
5-10-2011

QA SPECIALIST Cheeryl STRAUB AGENCY Tucson Police Department

DATE 05/05/2011 TIME 06:39

INTOXILYZER SERIAL # 80-001163 LOCATION TEAM 3

(☒) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(☒) 1. Clock time check.

(☒) 2. Date check.

OPERATIONAL TESTS

(☒) 1. Alcohol-free subject test result 0. 000 AC.

(☒) 2. Error recognition logic system functioning.

Not a Successfully Completed Test Sequence printed

(☒) 3. Proper sample recognition system.

Not a Successfully Completed Test Sequence printed

Deficient sample printed.

(☒) 4. Standard Calibration Check standard 0. 100 AC

Result: 0. 098 AC

Instrument is operating properly and accurately. YES ☒ NO ☐

COMMENTS

SIGNATURE

Cheeryl Straub

287 CS

B4

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001163
05/05/2011

06:40:33

Standard Lot#: 719919
Last Changed By: S. CIONE #50616

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SOAP, SOAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	06:41:27
Diagnostic Test	Pass	06:41:54
Air Blank	0.000	06:42:21
0.100 Cal Check	0.099	06:42:40
Air Blank	0.000	06:43:10
Subject Test	0.000	06:43:30
Air Blank	ABT*	06:43:37
Air Blank	0.000	06:44:05

*Sequence Aborted

Not a Successfully Completed Test Sequence

377 CS

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001163
05/05/2011

06:44:09

Bh

Standard Lot#: 719919
Last Changed By: S. CIONE #50616

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	06:44:41
Air Blank	IPS**	06:44:49

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

4g Fcs

154

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001163
05/05/2011

06:44:54

Standard Lot#: 719919
Last Changed By: S. CIONE #50616

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	06:45:42
Diagnostic Test	Pass	06:46:10
Air Blank	0.000	06:46:37
0.100 Cal Check	0.098	06:46:56
Air Blank	0.000	06:47:26
Subject Test	0.054*	06:50:34
Air Blank	0.000	06:51:07
Air Blank	0.000	06:51:34
Subject Test	XXX**	06:51:54
Air Blank	0.000	06:52:29

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

51705

154

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001163
05/05/2011 06:52:52

Standard Lot#: 719919
Last Changed By: S. CIONE #50616

QAS: C. STRAUB #50772
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	06:53:20
Cal Check	0.098	06:53:38
Air Blank	0.000	06:54:09



Operator's Signature

687 CS

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

TS4
5.10.2011

QA Specialist Cheryl STRAUB AGENCY Tucson Police Department
DATE 05/05/2011 TIME 06:39
Intoxilyzer Serial # 80-001163 Location TEAM 3

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution N/A AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Cheryl Straub

DPS Form Exh G-2 (Rev. 05-01)

727
CS


BS4

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001163
05/05/2011 06:39:02

Standard Lot#: 719919
Last Changed By: S. CIONE #50616

QAS: C. STRAUB #50772
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	06:39:29
Cal Check	0.099	06:39:48
Air Blank	0.000	06:40:18

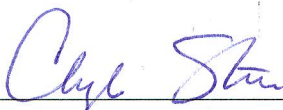

Operator's Signature