



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-001163 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4  
7-5-2011

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 06/30/11 Cal 1 0.097 AC Cal 2 0.097 AC

Between the dates of 06/03/11 and 06/30/11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.  
The instrument was returned to service on N/A.

Cheryl K. Straub

Cheryl K. Straub

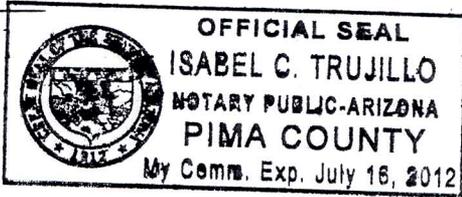
06/30/11

ADPS QAS Permit # 33408

SUBSCRIBED AND SWORN to before me this 30th day of June, 2011.

Isabel C. Trujillo

Notary Public



\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

177  
CS

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

TS4  
7-5-2011

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Cheryl STRAUB AGENCY Tucson Police Department

DATE 06/30/2011 TIME 06:35

INTOXILYZER SERIAL # 80-001163 LOCATION TEAM 3

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 097 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE Cheryl Straub

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

287  
CS

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-001163  
06/30/2011

06:36:46

TS4

Standard Lot#: 719919  
Last Changed By: S. CIONE #50616

QAS: C. STRAUB #50772  
TUCSON PD

Operator: C. STRAUB #50772  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	06:37:34
Diagnostic Test	Pass	06:38:02
Air Blank	0.000	06:38:29
0.100 Cal Check	0.097	06:38:48
Air Blank	0.000	06:39:18
Subject Test	0.000	06:39:35
Air Blank	ABT*	06:39:42
Air Blank	0.000	06:40:10

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

377  
CS

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-001163  
06/30/2011

06:40:15

154

Standard Lot#: 719919  
Last Changed By: S. CIONE #50616

QAS: C. STRAUB #50772  
TUCSON PD

Operator: C. STRAUB #50772  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	06:40:42
Air Blank	IPS**	06:40:50

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

497  
CS

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-001163  
06/30/2011

06:40:55

TS4

Standard Lot#: 719919  
Last Changed By: S. CIONE #50616

QAS: C. STRAUB #50772  
TUCSON PD

Operator: C. STRAUB #50772  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	06:41:41
Diagnostic Test	Pass	06:42:08
Air Blank	0.000	06:42:36
0.100 Cal Check	0.097	06:42:55
Air Blank	0.000	06:43:24
Subject Test	0.050*	06:46:33
Air Blank	0.000	06:47:05
Air Blank	0.000	06:47:32
Subject Test	XXX**	06:48:03
Air Blank	0.000	06:48:37

\*Deficient Sample  
\*\*Invalid Sample  
Begin new deprivation period

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

577  
CS

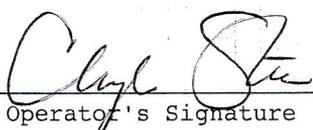
INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001163  
06/30/2011 06:49:17

B4

-----  
Standard Lot#: 719919  
Last Changed By: S. CIONE #50616  
-----

QAS: C. STRAUB #50772  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	06:49:45
Cal Check	0.097	06:50:03
Air Blank	0.000	06:50:33

  
-----  
Operator's Signature

697  
CS

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

TJG  
7.5.2011

QA Specialist Cheryl STRAUB AGENCY Tucson Police Department  
DATE 06/30/2011 TIME 06:35  
Intoxilyzer Serial # 80-001163 Location TEAM 3

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution N/A AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 097 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Cheryl Straub

DPS Form Exh G-2 (Rev. 05-01)

797  
08

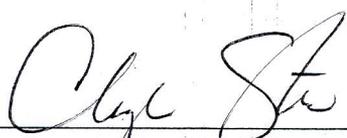
INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001163  
06/30/2011 06:35:19

B4

-----  
Standard Lot#: 719919  
Last Changed By: S. CIONE #50616  
-----

QAS: C. STRAUB #50772  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	06:35:47
Cal Check	0.097	06:36:05
Air Blank	0.000	06:36:35

  
-----  
Operator's Signature