



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-001163 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TSy  
8.2.2011

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 07/28/11 Cal 1 0.096 AC Cal 2 0.096 AC

Between the dates of 06/30/11 and 07/28/11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.  
The instrument was returned to service on N/A.

Cheryl K. Straub

Cheryl K. Straub


08/01/11

ADPS QAS Permit # 33408

SUBSCRIBED AND SWORN to before me this 1 day of August, 20 11.

Sue A. Soza

Notary Public

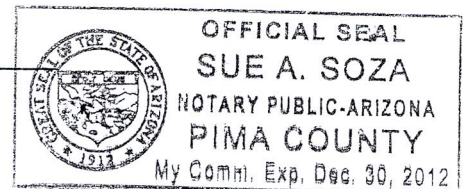

\*Includes documents specified in A.A.C. R13-10-104.A.5

188  
CS

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

TSg  
8.2.2011

QA SPECIALIST Cheeryl STRAUB AGENCY Tucson Police Department

DATE 07/28/2011 TIME 06:20

INTOXILYZER SERIAL # 80-001163 LOCATION TEAM 3

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.  
(✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.  
(✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed  
(✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.  
(✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 096 AC

Instrument is operating properly and accurately. YES ✓ NO       

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE

Cheeryl Straub

278  
CS

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-001163  
07/28/2011

06:22:17

154

Standard Lot#: 719919  
Last Changed By: S. CIONE #50616

QAS: C. STRAUB #50772  
TUCSON PD

Operator: C. STRAUB #50772  
TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	06:23:10
Diagnostic Test	Pass	06:23:37
Air Blank	0.000	06:24:05
0.100 Cal Check	0.096	06:24:23
Air Blank	0.000	06:24:53
Subject Test	0.000	06:25:07
Air Blank	ABT*	06:25:13
Air Blank	0.000	06:25:41

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

378cs

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-001163  
07/28/2011

06:25:46

BY

Standard Lot#: 719919  
Last Changed By: S. CIONE #50616

QAS: C. STRAUB #50772  
TUCSON PD

Operator: C. STRAUB #50772  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	06:26:13
Air Blank	IPS**	06:26:21

\*RFI Detect

\*\*Improper Sample

**Not a Successfully Completed Test Sequence**



498cs

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-001163  
07/28/2011

06:26:25

TS4

Standard Lot#: 719919  
Last Changed By: S. CIONE #50616

QAS: C. STRAUB #50772  
TUCSON PD

Operator: C. STRAUB #50772  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	06:27:10
Diagnostic Test	Pass	06:27:38
Air Blank	0.000	06:28:05
0.100 Cal Check	0.097	06:28:24
Air Blank	0.000	06:28:54
Subject Test	0.150*	06:32:02
Air Blank	0.000	06:32:34
Air Blank	0.000	06:33:01
Subject Test	0.099	06:33:28
Air Blank	0.000	06:34:01
Five Minute Wait	ABT**	
Air Blank	0.000	06:35:10

\*Deficient Sample

\*\*Sequence Aborted

→ Took sample as subject test and not mouth alcohol -  
Test Repeated - CS-50772 07/28/11

**Not a Successfully Completed Test Sequence**

588  
CS

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-001163  
07/28/2011

06:35:14

134

Standard Lot#: 719919  
Last Changed By: S. CIONE #50616

QAS: C. STRAUB #50772  
TUCSON PD

Operator: C. STRAUB #50772  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111

Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	06:36:04
Diagnostic Test	Pass	06:36:31
Air Blank	0.000	06:36:58
0.100 Cal Check	0.097	06:37:17
Air Blank	0.000	06:37:47
Subject Test	XXX*	06:38:07
Air Blank	0.000	06:38:41

\*Invalid Sample  
Begin new deprivation period

**Not a Successfully Completed Test Sequence**

688cs

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001163  
07/28/2011


06:39:09

154

-----  
Standard Lot#: 719919  
Last Changed By: S. CIONE #50616  
-----

QAS: C. STRAUB #50772  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	06:39:37
Cal Check	0.096	06:39:56
Air Blank	0.000	06:40:25

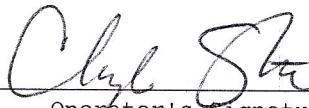

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Cheryl STRAUB AGENCY Tucson Police Department  
DATE 07/28/2011 TIME 06:20  
Intoxilyzer Serial # 80-001163 Location TEAM 3

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution N/A AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach  $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 096 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Cheryl Straub

DPS Form Exh G-2 (Rev. 05-01)



878cs

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001163  
07/28/2011


06:20:38

124

-----  
Standard Lot#: 719919  
Last Changed By: S. CIONE #50616  
-----

QAS: C. STRAUB #50772  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	06:21:06
Cal Check	0.096	06:21:24
Air Blank	0.000	06:21:54

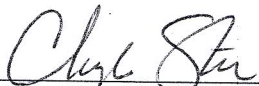

Operator's Signature