



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

754
11-7-2011

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001163 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 10-21-11 Cal 1 0.097 AC Cal 2 0.096 AC

Between the dates of 9-23-11 and 10-21-11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.


Clare Conley

Clare L. Conley

10-21-11

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 21, day of October, 20 11.

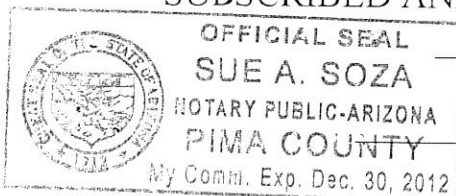

Sue A Soza

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Clare Conley AGENCY Tucson Police Department
DATE 10-21-11 TIME 0903
INTOXILYZER SERIAL # 80-001163 LOCATION ODM

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.
(✓) 2. Date check.

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0. 000 AC.
(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 096 AC

Instrument is operating properly and accurately. YES ✓ NO

COMMENTS None.

SIGNATURE Clare Conley

1567
On
TS4
11-7-2011

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001163
10/21/2011

09:04:25

287
Cu
TS4

Standard Lot#: 719919
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST,SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:05:14
Diagnostic Test	Pass	09:05:41
Air Blank	0.000	09:06:09
0.100 Cal Check	0.097	09:06:27
Air Blank	0.000	09:06:57
Subject Test	0.000	09:07:13
Air Blank	ABT*	09:07:19
Air Blank	RFI**	09:07:45

*Sequence Aborted

**RFI Detect

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001163
10/21/2011

09:07:50

Standard Lot#: 719919
Last Changed By: S. CIONE #50616
QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:08:37
Diagnostic Test	Pass	09:09:04
Air Blank	0.000	09:09:31
0.100 Cal Check	0.097	09:09:50
Air Blank	IPS*	09:09:58
Air Blank	0.000	09:10:26

*Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001163
10/21/2011

09:10:30

Standard Lot#: 719919
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:11:20
Diagnostic Test	Pass	09:11:47
Air Blank	0.000	09:12:15
0.100 Cal Check	0.097	09:12:33
Air Blank	0.000	09:13:03
Subject Test	0.176*	09:16:11
Air Blank	0.000	09:16:44
Air Blank	0.000	09:17:11
Subject Test	XXX**	09:17:37
Air Blank	0.000	09:18:11

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: ARIZONA

Serial Number: 80-001163

10/21/2011

09:18:33

Standard Lot#: 719919

Last Changed By: S. CIONE

#50616

QAS: C. CONLEY

#34517

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:19:01
0.100 Cal Check	0.096	09:19:20
Air Blank	0.000	09:19:50
End of test		

Clare Conley

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

6-267
CW

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TS4
11-7-2011

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Clare Conley AGENCY Tucson Police Department
DATE 10-21-11 TIME 0903
Intoxilyzer Serial # 80-001163 Location ODM

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 097 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001163
10/21/2011 09:03:04

Standard Lot#: 719919
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:03:31
Cal Check	0.097	09:03:50
Air Blank	0.000	09:04:20

Clare Conley
Operator's Signature