



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

SC
11/21/11

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001163 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 11/17/11 Cal 1 0.097 AC Cal 2 0.096 AC

Between the dates of 10/21/11 and 11/17/11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Cheryl K. Straub

Cheryl K. Straub


11/21/11

ADPS QAS Permit # 33408

SUBSCRIBED AND SWORN to before me this 21, day of November, 20 11.

A.M. Hernández

Notary Public

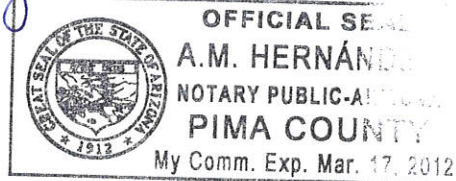

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

187
CS

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

SC
11/21/11

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Cheryl STRAUB AGENCY Tucson Police Department
DATE 11/17/11 TIME 06:36
INTOXILYZER SERIAL # 80-001163 LOCATION TEAM 3

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.
(✓) 2. Date check.

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0. 000 AC.
(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 096 AC

Instrument is operating properly and accurately. YES ☒ NO ☐

COMMENTS

SIGNATURE

Cheryl Straub

287cs

SC

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001163
11/17/2011

06:37:59

Standard Lot#: 719919
Last Changed By: S. CIONE #50616

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	06:38:50
Diagnostic Test	Pass	06:39:17
Air Blank	0.000	06:39:45
0.100 Cal Check	0.096	06:40:04
Air Blank	0.000	06:40:33
Subject Test	0.000	06:41:17
Air Blank	ABT*	06:41:24
Air Blank	0.000	06:41:52

*Sequence Aborted

Not a Successfully Completed Test Sequence

327
C8

SC

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001163
11/17/2011

06:41:59

Standard Lot#: 719919
Last Changed By: S. CIONE #50616

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	06:42:30
Air Blank	IPS**	06:42:38

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

487
CS

SC

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001163
11/17/2011

06:42:42

Standard Lot#: 719919
Last Changed By: S. CIONE #50616

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	06:43:29
Diagnostic Test	Pass	06:43:56
Air Blank	0.000	06:44:24
0.100 Cal Check	0.097	06:44:42
Air Blank	0.000	06:45:12
Subject Test	0.142*	06:48:21
Air Blank	0.000	06:48:53
Air Blank	0.000	06:49:20
Subject Test	XXX**	06:49:50
Air Blank	0.000	06:50:24

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

587
es


SC

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001163
11/17/2011 06:50:48

Standard Lot#: 719919
Last Changed By: S. CIONE #50616

QAS: C. STRAUB #50772
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	06:51:16
Cal Check	0.096	06:51:34
Air Blank	0.000	06:52:04

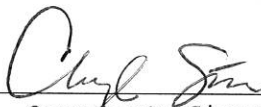

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

687
CS
SC
11/21/11

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Cheryl Straub AGENCY Tucson Police Department
DATE 11/17/11 TIME 06:36
Intoxilyzer Serial # 80-001163 Location TEAM 3

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution N/A AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 097 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Cheryl Straub

DPS Form Exh G-2 (Rev. 05-01)

777
CS

SC

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001163
11/17/2011 06:36:32

Standard Lot#: 719919
Last Changed By: S. CIONE #50616

QAS: C. STRAUB #50772
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	06:36:59
Cal Check	0.097	06:37:18
Air Blank	0.000	06:37:47



Operator's Signature