



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

TSK  
2.23.2012

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-001163 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 2/9/12 Cal 1 0.099 AC Cal 2 0.098 AC

Between the dates of 1/12/12 and 2/9/12 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


2/10/12

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 10<sup>th</sup>, day of February, 2012.

Elma Escarcega

Notary Public

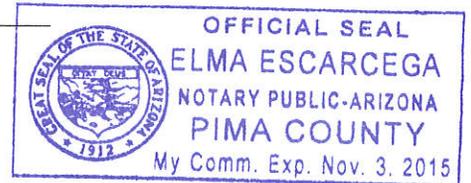

\*Includes documents specified in A.A.C. R13-10-104.A.5

1 of 7  
SC

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

TBY  
2.23.2012

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 2/9/12 TIME 1241

INTOXILYZER SERIAL # 80-001163 LOCATION ODM-Tm 3

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 098 AC

Instrument is operating properly and accurately. YES ✓ NO \_\_\_\_\_

COMMENTS

---



---



---



---

SIGNATURE Sarah Cione

2 of 7  
SC

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-001163  
02/09/2012 12:43:06

BM

---

Standard Lot#: 23511100A2  
Last Changed By: T. GALLEGOS #28847

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	12:43:54
Diagnostic Test	Pass	12:44:22
Air Blank	0.000	12:44:49
0.100 Cal Check	0.098	12:45:08
Air Blank	0.000	12:45:37
Subject Test	0.000	12:45:53
Air Blank	0.000	12:46:21
Five Minute Wait	ABT*	
Air Blank	0.000	12:47:11

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

3 of 7  
5C

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-001163  
02/09/2012 12:47:37

TS4

Standard Lot#: 23511100A2  
Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:48:22
Diagnostic Test	Pass	12:48:49
Air Blank	RFI*	12:48:59
Air Blank	IPS**	12:49:06

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-001163  
02/09/2012

12:49:16

4 of 7  
SC  
Tsh

---

Standard Lot#: 23511100A2  
Last Changed By: T. GALLEGOS #28847

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	12:50:14
Diagnostic Test	Pass	12:50:41
Air Blank	0.000	12:51:09
0.100 Cal Check	0.099	12:51:27
Air Blank	0.000	12:51:57
Subject Test	0.141*	12:55:05
Air Blank	0.000	12:55:37
Air Blank	0.000	12:56:05
Subject Test	XXX**	12:56:28
Air Blank	0.000	12:57:02

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

5 of 7 sc

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001163  
02/09/2012 12:57:16

TS4

-----  
Standard Lot#: 23511100A2  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	12:57:44
Cal Check	0.098	12:58:03
Air Blank	0.000	12:58:32

*Sarah Cione*

Operator's Signature

EXHIBIT G-2

6 of 7  
SC

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

TSC  
2.23.2012

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 2/9/12 TIME 1241  
Intoxilyzer Serial # 80-001163 Location ODM - Tr 3

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 099 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

7 of 7  
TSC

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001163  
02/09/2012 12:41:41

B4

-----  
Standard Lot#: 23511100A2  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	12:42:08
Cal Check	0.099	12:42:27
Air Blank	0.000	12:42:56

Sarah Cione  
Operator's Signature