



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

TS4
3-20-2012

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001163 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 3/9/12 Cal 1 0.100 AC Cal 2 0.098 AC

Between the dates of 2/9/12 and 3/9/12 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


3/9/12

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 9th, day of March, 2012.

Elma Escarcega

Notary Public

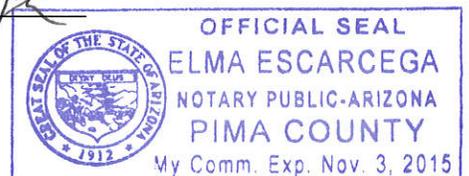

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

1 of 7
SC

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TSG
3-20-2012

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department
DATE 3/9/12 TIME 1045
INTOXILYZER SERIAL # 80-001163 LOCATION Tm 3-ODM

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 098 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

2 of 7
SC

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001163
03/09/2012

10:46:25

B4

Standard Lot#: 23511100A2
Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:47:09
Diagnostic Test	Pass	10:47:36
Air Blank	0.000	10:48:03
0.100 Cal Check	0.100	10:48:22
Air Blank	0.000	10:48:52
Subject Test	0.000	10:49:12
Air Blank	0.000	10:49:39
Five Minute Wait	ABT*	
Air Blank	0.000	10:50:16

*Sequence Aborted

Not a Successfully Completed Test Sequence

3 of 7
SC

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001163
03/09/2012

10:50:36

TS4

Standard Lot#: 23511100A2
Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	10:51:03
Air Blank	IPS**	10:51:11

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

4 of 7
SC

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001163
03/09/2012

10:51:21

TS4

Standard Lot#: 23511100A2
Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:52:04
Diagnostic Test	Pass	10:52:32
Air Blank	0.000	10:52:59
0.100 Cal Check	0.099	10:53:18
Air Blank	0.000	10:53:47
Subject Test	0.042*	10:56:56
Air Blank	0.000	10:57:28
Air Blank	0.000	10:57:55
Subject Test	XXX**	10:58:26
Air Blank	0.000	10:58:59

*Deficient Sample
 **Invalid Sample
 Begin new deprivation period

Not a Successfully Completed Test Sequence

5 of 7
SC

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001163
03/09/2012 10:59:16

TS4

Standard Lot#: 23511100A2
Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:59:44
Cal Check	0.098	11:00:02
Air Blank	0.000	11:00:32

Sarah Cione

Operator's Signature

6 of 7
5c

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TS4
3.20.2012

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 3/9/12 TIME 1045
Intoxilyzer Serial # 80-001163 Location Tm3-ODM

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001163
03/09/2012 10:45:00

7 of 7
50

Standard Lot#: 23511100A2
Last Changed By: T. GALLEGOS #28847

BSH

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:45:27
Cal Check	0.100	10:45:46
Air Blank	0.000	10:46:15

Sarah Cione

Operator's Signature