



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

SC
6/18/12

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001163 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 5.31.2012 Cal 1 0.098 AC Cal 2 0.098 AC

Between the dates of 5.3.2012 and 5.31.2012 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.
The instrument was returned to service on —.

Terry S. Gallegos

Terry S. Gallegos


6.13.2012

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 13th, day of June, 2012.

Elma Escarcega

Notary Public

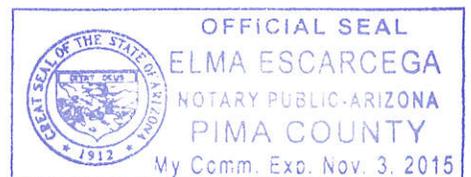

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

10 of 7
TS4
SC
6/18/12

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 5.31.2012 TIME 1041
INTOXILYZER SERIAL # 80-001163 LOCATION ODM

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 098 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Terry Gallegos

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001163
05/31/2012

10:43:16

2057
TS4

Standard Lot#: 23511100A2
Last Changed By: T. GALLEGOS #28847

SC

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:44:09
Diagnostic Test	Pass	10:44:37
Air Blank	0.000	10:45:04
0.100 Cal Check	0.099	10:45:23
Air Blank	0.000	10:45:52
Subject Test	0.000	10:46:34
Air Blank	0.000	10:47:02
Five Minute Wait	ABT*	
Air Blank	0.000	10:47:35

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001163
05/31/2012 10:48:00

3057
TS4
SC

Standard Lot#: 23511100A2
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:48:58
Diagnostic Test	Pass	10:49:25
Air Blank	0.000	10:49:53
0.100 Cal Check	RFI*	10:50:08
Air Blank	IPS**	10:50:20

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-001163
05/31/2012

10:50:36

Standard Lot#: 23511100A2
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST,SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:51:27
Diagnostic Test	Pass	10:51:55
Air Blank	0.000	10:52:22
0.100 Cal Check	0.098	10:52:41
Air Blank	0.000	10:53:10
Subject Test	0.172*	10:56:19
Air Blank	0.000	10:56:51
Air Blank	0.000	10:57:18
Subject Test	XXX**	10:57:42
Air Blank	0.000	10:58:15

*Deficient Sample

**Invalid Sample

Begin new deprivation period

4057
TD4
SC

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001163
05/31/2012 10:59:34

Standard Lot#: 23511100A2
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:00:01
Cal Check	0.098	11:00:20
Air Blank	0.000	11:00:50

5087
TB4
SC



Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

6077
TS4
SC
6/18/12

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 5.31.2012 TIME 1041
Intoxilyzer Serial # 80-001163 Location ODM

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 098 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)


INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001163
05/31/2012 10:41:44

Standard Lot#: 23511100A2
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Fof7
B4
SC

Test	g/210L	Time
Air Blank	0.000	10:42:11
Cal Check	0.098	10:42:30
Air Blank	0.000	10:42:59

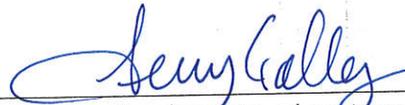

Operator's Signature