



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

SC  
4/19/13

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-001163 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 4-3-2013 Cal 1 0.097 AC Cal 2 0.097 AC

Between the dates of 3-7-2013 and 4-3-2013 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.  
The instrument was returned to service on —.


Terry S. Gallegos

4-16-2013

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 16<sup>th</sup>, day of April, 2013.

Notary Public

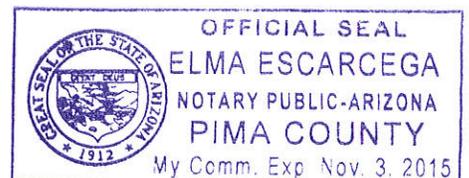

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7  
DH  
SC  
4/19/13

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 4.3.2013 TIME 1724

INTOXILYZER SERIAL # 80-001163 LOCATION CrimeLab

(  ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (  ) 1. Clock time check.
- (  ) 2. Date check.

OPERATIONAL TESTS

- (  ) 1. Alcohol-free subject test result 0.000 AC.
- (  ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (  ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (  ) 4. Standard Calibration Check standard 0.100 AC  
Result: 0.097 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE Terry Gallegos

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001163  
Core Version: 8105.48  
04/03/2013

17:28:54

2077  
DJ4

5<

---

Standard Lot#: 23511100A2  
Last Changed By: T. GALLEGOS #28847

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	17:29:55
Diagnostic Test	Pass	17:30:30
Air Blank	0.000	17:30:57
0.100 Cal Check	0.098	17:31:18
Air Blank	0.000	17:31:48
Subject Test	0.000	17:32:08
Air Blank	0.000	17:32:36
Wait	ABT*	
Air Blank	0.000	17:33:10

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001163  
Core Version: 8105.48  
04/03/2013

17:33:28

30f7  
TD4

5C

---

Standard Lot#: 23511100A2  
Last Changed By: T. GALLEGOS #28847

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	RFI*	17:34:06
Air Blank	IPS**	17:34:14

\*RFI Detect

\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001163  
Core Version: 8105.48  
04/03/2013

17:34:21

40f7  
TD4  
SC

---

Standard Lot#: 23511100A2  
Last Changed By: T. GALLEGOS #28847

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST, SQAP

DOB: 11/11/1987

Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	17:35:11
Diagnostic Test	Pass	17:35:48
Air Blank	0.000	17:36:15
0.100 Cal Check	0.097	17:36:36
Air Blank	0.000	17:37:05
Subject Test	0.316*	17:40:13
Air Blank	0.000	17:40:47
Air Blank	0.000	17:41:14
Subject Test	XXX**	17:42:00
Air Blank	0.000	17:42:34

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001163  
Core Version: 8105.48  
04/03/2013 18:12:07

50f7  
D4

-----  
Standard Lot#: 23511100A2  
Last Changed By: T. GALLEGOS #28847  
-----

SC

-----  
QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

-----  
Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	18:12:37
0.100 Cal Check	0.097	18:12:58
Air Blank	0.000	18:13:28

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

6077  
TJG

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

SC  
4/19/13

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 2887 AGENCY Tucson Police Department  
DATE 4.3.2013 TIME 1724  
Intoxilyzer Serial # 80-001163 Location CrimeLab

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 097 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001163  
Core Version: 8105.48  
04/03/2013 17:24:17

7 of 7  
TD4

-----  
Standard Lot#: 23511100A2  
Last Changed By: T. GALLEGOS #28847  
-----

SC

-----  
QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

-----  
Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	17:24:48
0.100 Cal Check	0.097	17:25:08
Air Blank	0.000	17:25:38
Air Blank	0.000	17:28:06

\*Input Timed Out