



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

SC
6/6/13

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001163 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 5.30.2013 Cal 1 0.097 AC Cal 2 0.095 AC

Between the dates of 5.1.2013 and 5.30.2013 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.
The instrument was returned to service on —.

Terry S. Gallegos

Terry S. Gallegos


6-4-2013

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 2nd, day of June, 2013.

Elma Escarcega

Notary Public

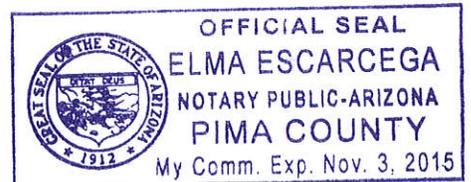

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7
TS4

SC
6/6/13

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 5.30.2013 TIME 0939

INTOXILYZER SERIAL # 80-001163 LOCATION ODE

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 095 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE 

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001163
Core Version: 8105.48
05/30/2013

09:41:01

2077
TS4
SC

Standard Lot#: 23511100A2
Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:41:52
Diagnostic Test	Pass	09:42:27
Air Blank	0.000	09:42:54
0.100 Cal Check	0.096	09:43:15
Air Blank	0.000	09:43:45
Subject Test	0.000	09:44:12
Air Blank	0.000	09:44:40
Wait	ABT*	
Air Blank	0.000	09:45:21

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001163
Core Version: 8105.48
05/30/2013 09:45:39

3077
TS4
5C

Standard Lot#: 23511100A2
Last Changed By: T. GALLEGOS #28847

QAS: S. CIONE #50616
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:46:31
Diagnostic Test	Pass	09:47:06
Air Blank	0.000	09:47:34
0.100 Cal Check	RFI*	09:47:49
Air Blank	IPS**	09:47:59

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001163
Core Version: 8105.48
05/30/2013 09:48:04

4 of 7
TS4

Standard Lot#: 23511100A2
Last Changed By: T. GALLEGOS #28847

SC

QAS: S. CIONE #50616
TUCSON PD

Operator: T. GALLEGOS #28847
TPD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:49:18
Diagnostic Test	Pass	09:49:53
Air Blank	0.000	09:50:21
0.100 Cal Check	0.096	09:50:41
Air Blank	0.000	09:51:11
Subject Test	0.096*	09:54:18
Air Blank	0.000	09:54:50
Air Blank	0.000	09:55:18
Subject Test	XXX**	09:55:39
Air Blank	0.000	09:56:13

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001163
Core Version: 8105.48
05/30/2013 09:56:42

Soft
TSh
SC

Standard Lot#: 23511100A2
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:57:13
0.100 Cal Check	0.095	09:57:34
Air Blank	0.000	09:58:03

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

6077
TS4

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

SC
6/6/13

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 5.30.2013 TIME 0939
Intoxilyzer Serial # 80-001163 Location ODE

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 097 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001163
Core Version: 8105.48
05/30/2013 09:39:12

Fof7
TS4

Standard Lot#: 23511100A2
Last Changed By: T. GALLEGOS #28847

SC

QAS: S. CIONE #50616
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	09:39:39
0.100 Cal Check	0.097	09:40:00
Air Blank	0.000	09:40:29