



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

SC  
8/8/13

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-001163 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 7.25.2013 Cal 1 0.096 AC Cal 2 0.095 AC

Between the dates of 6.27.2013 and 7.25.2013 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on         .  
The instrument was returned to service on         .

Terry Gallegos

Terry S. Gallegos


8.7.2013

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 7<sup>th</sup>, day of August, 2013.

Elma Escarcega

Notary Public

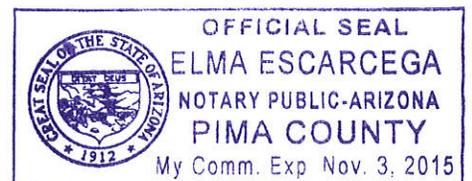

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7

TSh

SC  
8/8/13

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 7.25.2013 TIME 0915

INTOXILYZER SERIAL # 80001163 LOCATION ODE

(  ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (  ) 1. Clock time check.
- (  ) 2. Date check.

OPERATIONAL TESTS

- (  ) 1. Alcohol-free subject test result 0. 000 AC.
- (  ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (  ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (  ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 095 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE Terry Gallegos

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001163  
Core Version: 8105.48  
07/25/2013

09:17:12

2087  
TD4

Standard Lot#: 23511100A2  
Last Changed By: T. GALLEGOS #28847

5<

OAS: T. GALLEGOS #28847  
TUCSON PD

Operator: T. GALLEGOS #28847  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:18:07
Diagnostic Test	Pass	09:18:42
Air Blank	0.000	09:19:09
0.100 Cal Check	0.096	09:19:30
Air Blank	0.000	09:19:59
Subject Test	0.000	09:20:38
Air Blank	0.000	09:21:05
Wait	ABT*	
Air Blank	0.000	09:21:39

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001163  
Core Version: 8105.48  
07/25/2013

09:22:09

30 of 7  
T24  
5C

---

Standard Lot#: 23511100A2  
Last Changed By: T. GALLEGOS #28847

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST,SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	09:23:06
Diagnostic Test	Pass	09:23:41
Air Blank	0.000	09:24:08
0.100 Cal Check	RFI*	09:24:21
Air Blank	IPS**	09:24:30

\*RFI Detect

\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001163  
Core Version: 8105.48  
07/25/2013

09:24:45

40f7  
TDG

5C

---

Standard Lot#: 23511100A2  
Last Changed By: T. GALLEGOS #28847

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST,SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	09:25:40
Diagnostic Test	Pass	09:26:15
Air Blank	0.000	09:26:42
0.100 Cal Check	0.095	09:27:03
Air Blank	0.000	09:27:33
Subject Test	0.104*	09:30:40
Air Blank	0.000	09:31:12
Air Blank	0.000	09:31:39
Subject Test	XXX**	09:32:22
Air Blank	0.000	09:32:56

\*Deficient Sample  
\*\*Invalid Sample  
Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001163  
Core Version: 8105.48  
07/25/2013 09:33:16

Sof 7  
TX

-----  
Standard Lot#: 23511100A2  
Last Changed By: T. GALLEGOS #28847

SC

-----  
QAS: T. GALLEGOS #28847  
TUCSON PD

-----  
Operator: T. GALLEGOS #28847  
TUCSON PD

-----  
Test g/210L Time  
-----  
Air Blank 0.000 09:33:47  
0.100 Cal Check 0.095 09:34:08  
Air Blank 0.000 09:34:37

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

60f7  
TS4  
SC  
8/8/13

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 7.25.2013 TIME 0915  
Intoxilyzer Serial # 80-001163 Location ODE

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 096 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001163  
Core Version: 8105.48  
07/25/2013 09:15:14

Fof 7  
TD4

8/8/13 SC

-----  
Standard Lot#: 23511100A2  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:15:55
0.100 Cal Check	0.096	09:16:15
Air Blank	0.000	09:16:45

Preventative Maintenance NOT  
Performed -- Timer NOT Reset