



CERTIFICATION by CUSTODIAN of RECORDS
 Tucson Police Department Crime Laboratory
 1306 W. Miracle Mile, Tucson, AZ 85705

SC
 12/11/13

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001163 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 11.14.2013 Cal 1 0.097 AC Cal 2 0.096 AC

Between the dates of 10.17.2013 and 11.14.2013 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.
 The instrument was returned to service on —.

Terry Gallegos

Terry S. Gallegos


12.3.2013

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 3rd day of December, 2013.

Elma Escarcega

Notary Public

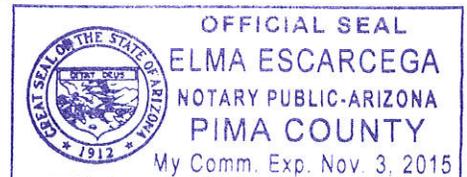

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

10/10
TDH
SC
12/11/13

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 11-14-2013 TIME 1703

INTOXILYZER SERIAL # 80-001163 LOCATION Crime Lab

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 096 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE 

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001163
Core Version: 8105.48
11/14/2013

17:12:56

20/10
TS
SC

Standard Lot#: 30412100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	17:14:00
Diagnostic Test	Pass	17:14:35
Air Blank	0.000	17:15:02
0.100 Cal Check	0.096	17:15:23
Air Blank	0.000	17:15:52
Subject Test	0.000	17:16:29
Air Blank	0.000	17:16:56
Wait	ABT*	
Air Blank	0.000	17:17:31

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001163
Core Version: 8105.48
11/14/2013

17:18:03

30f10
TS4

5C

Standard Lot#: 30412100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	17:18:51
Diagnostic Test	Pass	17:19:26
Air Blank	0.000	17:19:54
0.100 Cal Check	RFI*	17:20:09
Air Blank	IPS**	17:20:16

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001163
Core Version: 8105.48
11/14/2013

17:30:29

40 of 10
B4
5C

Standard Lot#: 30412100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	17:31:27
Diagnostic Test	Pass	17:32:02
Air Blank	0.000	17:32:30
0.100 Cal Check	0.096	17:32:50
Air Blank	0.000	17:33:20
Subject Test	0.016*	17:36:27
Air Blank	0.000	17:37:00
Air Blank	0.000	17:37:27
Reference	RFI**	17:37:34
Air Blank	0.000	17:38:03

*Deficient Sample

**RFI Detect

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001163
Core Version: 8105.48
11/14/2013

17:39:03

50f10
T36

5c

Standard Lot#: 30412100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	17:40:05
Diagnostic Test	Pass	17:40:40
Air Blank	0.000	17:41:07
0.100 Cal Check	0.096	17:41:28
Air Blank	0.000	17:41:58
Subject Test	XXX*	17:42:19
Air Blank	0.000	17:42:53

*Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001163
Core Version: 8105.48
11/14/2013 17:43:24

60f 10
TS4

Standard Lot#: 30412100A1
Last Changed By: S. CIONE #50616

SC

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	17:43:53
0.100 Cal Check	0.096	17:44:14
Air Blank	0.000	17:44:43

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

Fof 10
TS4

SC
12/11/13

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 11-14-2013 TIME 1703
Intoxilyzer Serial # 80-001163 Location CrimeLab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 097 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001163
Core Version: 8105.48
11/14/2013 17:03:19

80f10
TS4

Standard Lot#: 30412100A1
Last Changed By: S. CIONE #50616

5c

QAS: S. CIONE #50616
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	17:03:46
0.100 Cal Check	0.097	17:04:07
Air Blank	0.000	17:04:36



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

9 of 10
Tdy
sc

INTOXILYZER SERVICE RECORD
FOR
CALIBRATION GAS TANK CHANGE

Date: 10/31/13
QAS: Sarah Cione
Intoxilyzer Serial Number: 80-0G1163
New Lot Number: 30412100A1
New Expiration Date: 12/1/2014

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."
Select "D" for dry gas and enter the following information:
 - Enter Target Value: 0.100
 - Enter Control Lot #: From the tank label
 - Enter Exp Date: From the tank label (mm/dd/yyyy)
 - Edit info Y/N? Y or N
 - Update Std Info? Y ENTER
 - "Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing."
Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 098 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-001163
Core Version: 8105.48
10/31/2013 15:01:54

10 of 10
TSY
SC

Standard Lot#: 30412100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	15:02:22
0.100 Cal Check	0.098	15:02:43
Air Blank	0.000	15:03:13

Preventative Maintenance
Performed -- Timer Reset