



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

SC  
1/14/14

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-001163 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 1-9-2014 Cal 1 0.096 AC Cal 2 0.096 AC

Between the dates of 12-13-2013 and 1-9-2014 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.  
The instrument was returned to service on —.

Terry S. Gallegos

Terry S. Gallegos


1-13-2014

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 13<sup>th</sup>, day of January, 2014.

Elma Escarcega

Notary Public

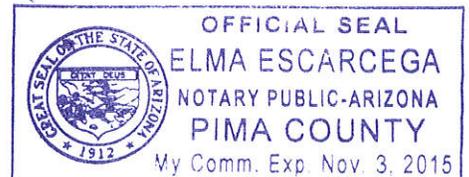

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

1 of 7  
TSg  
SC  
1/14/14

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 1-9-2014 TIME 1733  
INTOXILYZER SERIAL # 80-001163 LOCATION CrimeLab

(  ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (  ) 1. Clock time check.
- (  ) 2. Date check.

OPERATIONAL TESTS

- (  ) 1. Alcohol-free subject test result 0. 000 AC.
- (  ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (  ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (  ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 096 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE

Terry Gallegos

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001163  
Core Version: 8105.48  
01/09/2014

17:39:42

20f7  
TSh  
5C

---

Standard Lot#: 30412100A1  
Last Changed By: S. CIONE #50616

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	17:40:38
Diagnostic Test	Pass	17:41:13
Air Blank	0.000	17:41:40
0.100 Cal Check	0.096	17:42:01
Air Blank	0.000	17:42:31
Subject Test	0.000	17:43:46
Air Blank	0.000	17:44:14
Wait	ABT*	
Air Blank	0.000	17:44:48

\*Sequence Aborted

### Not a Successfully Completed Test Sequence

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001163  
Core Version: 8105.48  
01/09/2014

17:45:14

30f7  
734  
56

---

Standard Lot#: 30412100A1  
Last Changed By: S. CIONE #50616

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	17:46:36
Diagnostic Test	Pass	17:47:11
Air Blank	0.000	17:47:38
0.100 Cal Check	RFI*	17:47:52
Air Blank	IPS**	17:48:00

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001163  
Core Version: 8105.48  
01/09/2014

17:48:52

40f7  
TS4

5C

---

Standard Lot#: 30412100A1  
Last Changed By: S. CIONE #50616

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST,SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	17:49:39
Diagnostic Test	Pass	17:50:14
Air Blank	0.000	17:50:41
0.100 Cal Check	0.096	17:51:02
Air Blank	0.000	17:51:32
Subject Test	0.032*	17:54:39
Air Blank	0.000	17:55:11
Air Blank	0.000	17:55:38
Subject Test	XXX**	17:55:55
Air Blank	0.000	17:56:28

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001163  
Core Version: 8105.48  
01/09/2014 17:57:28

50f7  
TS4

-----  
Standard Lot#: 30412100A1  
Last Changed By: S. CIONE #50616  
-----

5C

-----  
QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

-----  
Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	17:57:57
0.100 Cal Check	0.096	17:58:18
Air Blank	0.000	17:58:47

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

60f7  
TSh  
SC  
1/14/14

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 1-9-2014 TIME 1733  
Intoxilyzer Serial # 80-001163 Location CrimeLab

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 096 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE 

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-001163  
Core Version: 8105.48  
01/09/2014 17:33:59

70f7  
TS4  
5C

-----  
Standard Lot#: 30412100A1  
Last Changed By: S. CIONE #50616  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD CRIME LAB  
-----

Test	g/210L	Time
Air Blank	0.000	17:34:27
0.100 Cal Check	0.096	17:34:48
Air Blank	0.000	17:35:17