



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

SC  
7/8/14

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-001163 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6.26.2014 Cal 1 0.097 AC Cal 2 0.096 AC

Between the dates of 5.29.2014 and 6.26.2014 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
The instrument was returned to service on -.

Terry S. Gallegos

Terry S. Gallegos


7-2-2014

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 2nd, day of July, 2014.

Elma Escarcega

Notary Public

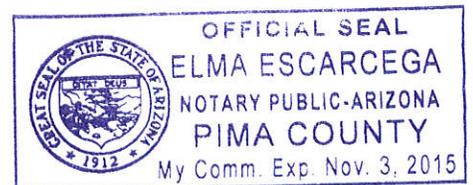

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

10f7  
TSy

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

sc  
7/8/14

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 6.26.2014 TIME 1200

INTOXILYZER SERIAL # 80-001163 LOCATION Sgt's Office ODM

(  ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (  ) 1. Clock time check.
- (  ) 2. Date check.

OPERATIONAL TESTS

- (  ) 1. Alcohol-free subject test result 0. 000 AC.
- (  ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (  ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (  ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 096 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS \_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Terry Gallegos

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

2017  
TS4

50

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-001163  
Core Version: 8105.48  
06/26/2014 12:02:40

-----  
Standard Lot#: 25913100A2  
Last Changed By: S. CIONE #50616  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Subject: TEST, SOAP  
DOB: 11/11/1987  
Sex: F Weight: 111  
-----

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:03:33
Diagnostic Test	Pass	12:04:10
Air Blank	0.000	12:04:38
0.100 Cal Check	0.097	12:04:58
Air Blank	0.000	12:05:28
Subject Test	0.000	12:06:06
Air Blank	0.000	12:06:33
Wait	ABT*	
Air Blank	0.000	12:07:12

\*Sequence Aborted

Not a Successfully  
Completed  
Test Sequence

|

3 of 7  
DJ  
SC

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-001163  
Core Version: 8105.48  
06/26/2014 12:08:24

-----  
Standard Lot#: 25913100A2  
Last Changed By: S. CIONE #50616  
-----

OAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Subject: TEST,SOAP  
DOB: 11/11/1987  
Sex: F Weight: 111  
-----

15 Minute Deprivation Period? Yes  
-----

Test	g/210L	Time
Air Blank	0.000	12:09:13
Diagnostic Test	Pass	12:09:51
Air Blank	0.000	12:10:18
0.100 Cal Check	RFI*	12:10:32
Air Blank	IPS**	12:10:40

\*RFI Detect  
\*\*Improper Sample

Not a Successfully  
Completed  
Test Sequence



4 of 7  
DJ  
SC

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-001163  
Core Version: 8105.48  
06/26/2014 12:10:50

-----  
Standard Lot#: 25913100A2  
Last Changed By: S. CIONE #50616  
-----

DAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Subject: TEST,SOAP  
DOB: 11/11/1987  
Sex: F Weight: 111  
-----

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:11:45
Diagnostic Test	Pass	12:12:22
Air Blank	0.000	12:12:50
0.100 Cal Check	0.096	12:13:10
Air Blank	0.000	12:13:40
Subject Test	0.365*	12:16:47
Air Blank	0.000	12:17:20
Air Blank	0.000	12:17:47
Subject Test	XXX**	12:18:15
Air Blank	0.000	12:18:48

\*Deficient Sample  
\*\*Invalid Sample  
Begin new deprivation period

Not a Successfully  
Completed  
Test Sequence

Sof 7  
TD4

SC

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-001163  
Core Version: 8105.48  
06/26/2014 12:19:14

-----  
Standard Lot#: 25913100A2  
Last Changed By: S. CIONE #50616  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	12:19:43
0.100 Cal Check	0.096	12:20:04
Air Blank	0.000	12:20:33

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6077  
TS4  
5C  
7/8/14

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 6-26-2014 TIME 1200  
Intoxilyzer Serial # 80-001163 Location Sgt's office ODM

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C + 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 097 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Instrument Initialization  
11:19:31 06/26/2014

7 of 7  
TS4  
58

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-001163  
Core Version: 8105.48  
06/26/2014 12:00:48

-----  
Standard Lot#: 25913100A2  
Last Changed By: S. CIONE #50616  
-----

OAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD CRIME LAB  
-----

Test	g/210L	Time
Air Blank	0.000	12:01:15
0.100 Cal Check	0.097	12:01:36
Air Blank	0.000	12:02:06