



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

SC

10/13/14

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001163 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 9.18.2014 Cal 1 0.095 AC Cal 2 0.094 AC

Between the dates of 8.21.2014 and 9.18.2014 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.
The instrument was returned to service on —.


Terry S. Gallegos

Terry S. Gallegos

10.9.2014

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 9th day of October, 2014.

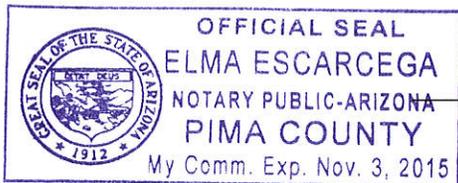

Elma Escarcega

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 9
B4

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

SC
10/13/14

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 9-18-2014 TIME 1656

INTOXILYZER SERIAL # 80-001163 LOCATION ODM Sgt's office

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 094 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE 

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

205 9
754
SC

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-001163
Core Version: 8105.48
09/18/2014 16:58:45

Standard Lot#: 25913100A2
Last Changed By: S. CIONE #50616

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:59:42
Diagnostic Test	Pass	17:00:19
Air Blank	0.000	17:00:47
0.100 Cal Check	0.095	17:01:08
Air Blank	0.000	17:01:37
Subject Test	0.000	17:02:23
Air Blank	0.000	17:02:51
Wait	ABT*	
Air Blank	0.000	17:03:31

*Sequence Aborted

Not a Successfully
Completed
Test Sequence

3089
784
SC

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-001163
Core Version: 8105.48
09/18/2014 17:05:12

Standard Lot#: 25913100A2
Last Changed By: S. CIONE #50616

DAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST,SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	17:06:50
Diagnostic Test	Pass	17:07:27
Air Blank	0.000	17:07:55
0.100 Cal Check	RFI*	17:08:08
Air Blank	IPS**	17:08:17

*RFI Detect
**Improper Sample

Not a Successfully
Completed
Test Sequence

4 of 9
TS4
SC

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-001163
Core Version: 8105.48
09/18/2014 17:08:48

Standard Lot#: 25913100A2
Last Changed By: S. CIONE #50616

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SOAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	17:09:38
Diagnostic Test	Pass	17:10:15
Air Blank	0.000	17:10:42
0.100 Cal Check	0.095	17:11:03
Air Blank	0.000	17:11:33
Subject Test	0.204*	17:14:40
Air Blank	0.000	17:15:13
Air Blank	0.000	17:15:40
Subject Test	XXX**	17:16:07
Air Blank	0.000	17:16:41

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully
Completed
Test Sequence

5089
154

SC

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-D01163
Core Version: 8105.48
09/18/2014 17:17:12

Standard Lot#: 25913100A2
Last Changed By: S. CIONE #50616

DAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	17:17:41
0.100 Cal Check	0.094	17:18:02
Air Blank	0.000	17:18:32

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

6059
B4
SC
10/13/14

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 9-18-2014 TIME 1656
Intoxilyzer Serial # 80-001163 Location ODM Sgt's Office

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 095 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE 

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Instrument Initialization
16:19:03 09/18/2014

Fof 9
Tdy
sc

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-001163
Core Version: 8105.48
09/18/2014 16:56:56

Standard Lot#: 25913100A2
Last Changed By: S. CIONE #50616

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	16:57:25
0.100 Cal Check	0.095	16:57:46
Air Blank	0.000	16:58:15

Preventative Maintenance
Performed -- Timer Reset



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

8059
TS4

INTOXILYZER SERVICE RECORD
FOR
CALIBRATION GAS TANK CHANGE

SC
10/3/14

Date: 9.18.2014
QAS: Terry Gallegos
Intoxilyzer Serial Number: 80-001163
New Lot Number: 115141002A
New Expiration Date: 05.01.2016

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."
Select "D" for dry gas and enter the following information:
 - Enter Target Value: 0.100
 - Enter Control Lot #: From the tank label
 - Enter Exp Date: From the tank label (mm/dd/yyyy)
 - Edit info Y/N? Y or N
 - Update Std Info? Y ENTER
 - "Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing."
Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 098 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

9 of 9
TS6

56

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-001163
Core Version: 8105.48
09/18/2014 17:37:35

Standard Lot#: 115141002A
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	17:38:05
0.100 Cal Check	0.098	17:38:26
Air Blank	0.000	17:38:55

Preventative Maintenance
Performed -- Timer Reset