



CERTIFICATION by CUSTODIAN of RECORDS
 Tucson Police Department Crime Laboratory
 1306 W. Miracle Mile, Tucson, AZ 85705

PH
 4.13.2015

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001163 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 4/1/2015 Cal 1 0.098 AC Cal 2 0.097 AC

Between the dates of 3/5/2015 and 4/1/2015 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on _____.
 The instrument was returned to service on _____.

William L Loomis

William L. Loomis


4/2/2015

ADPS QAS Permit # 33720

SUBSCRIBED AND SWORN to before me this 2nd, day of April, 2015.

Edna Escarcega

Notary Public

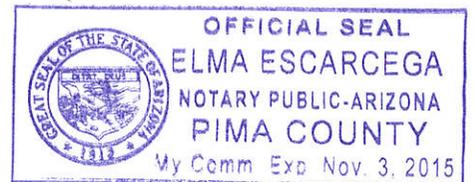

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 8
WLL
By
4-13-2015

QA SPECIALIST WILL LOOMIS 53545 AGENCY Tucson Police Department

DATE 4/1/2015 TIME 17:16

INTOXILYZER SERIAL # 80-001163 LOCATION CRIME LAB

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 097 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Will Loomis

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001163
Core Version: 8105.48
04/01/2015 17:19:34

2 of 8
WLV
B4

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SQAP,TEST
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	17:20:44
Diagnostic Test	Pass	17:21:19
Air Blank	0.000	17:21:46
0.100 Cal Check	0.097	17:22:07
Air Blank	0.000	17:22:37
Subject Test	0.000	17:23:38
Air Blank	0.000	17:24:06
Wait	ABT*	
Air Blank	0.000	17:26:00

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001163
Core Version: 8105.48
04/01/2015 17:26:16

3 of 8
Wkt
TS4

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SQAP,TEST
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	17:27:07
Diagnostic Test	Pass	17:27:42
Air Blank	0.000	17:28:10
0.100 Cal Check	RFI*	17:28:25
Air Blank	IPS**	17:28:34

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001163
Core Version: 8105.48
04/01/2015 17:28:40

4 of 8
WLL
TBH

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	17:29:33
Diagnostic Test	Pass	17:30:08
Air Blank	0.000	17:30:35
0.100 Cal Check	0.097	17:30:56
Air Blank	0.000	17:31:26
Subject Test	RGE*	17:31:44
Air Blank	0.000	17:32:20

*Range Exceeded

ON DEFICIENT SAMPLE ATTEMPT,

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001163
Core Version: 8105.48
04/01/2015 17:32:24

5 of 8
WLC
134

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SQAP,TEST
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	17:33:14
Diagnostic Test	Pass	17:33:50
Air Blank	0.000	17:34:17
0.100 Cal Check	0.098	17:34:38
Air Blank	0.000	17:35:07
Subject Test	0.327*	17:38:15
Air Blank	0.000	17:38:47
Air Blank	0.000	17:39:15
Subject Test	XXX**	17:39:31
Air Blank	0.000	17:40:04

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD CALIBRATION CHECK PROCEDURE

7 of 8
TS4
WLL
4.13.2015

QA Specialist Will Loomis 53545 AGENCY Tucson Police Department
DATE 4/1/2015 TIME 17:16
Intoxilyzer Serial # 80-001163 Location CRIME LAB

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 098 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Will Loomis

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-001163

Core Version: 8105.48

04/01/2015

17:16:53

8 of 8
WLL
TBH

Standard Lot#: 11514100A2

Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545

TUCSON PD

Operator: W. LOOMIS #53545

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	17:17:56
0.100 Cal Check	0.098	17:18:17
Air Blank	0.000	17:18:47

Preventative Maintenance
Performed -- Timer Reset