



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

TS4  
6-1-2015

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-001163 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 5-29-15 Cal 1 0.096 AC Cal 2 0.096 AC

Between the dates of 4-30-15 and 5-29-15 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.  
The instrument was returned to service on —.


Amy L. Autret

6-1-15

ADPS QAS Permit # 33718

SUBSCRIBED AND SWORN to before me this 1st, day of June, 2015.

Notary Public

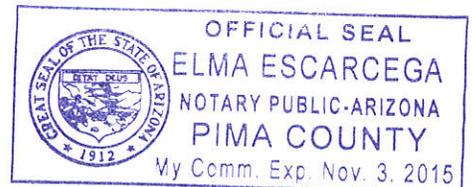

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

1018  
JL  
TSY  
6.1.2015

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Amy Autret 53546 AGENCY Tucson Police Department

DATE 5-29-15 TIME 0741

INTOXILYZER SERIAL # 80-001163 LOCATION Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS


- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 096 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE

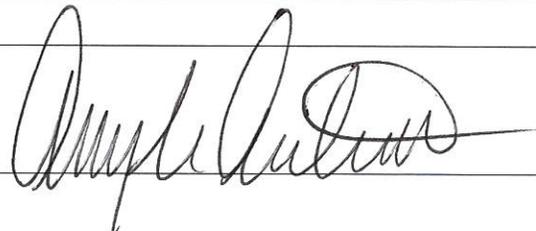

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001163  
Core Version: 8105.48  
05/29/2015 07:44:41

2068  
JL  
TS4

---

Standard Lot#: 11514100A2  
Last Changed By: T. GALLEGOS #28847

---

QAS: A. AUTRET #53546  
TUCSON PD

---

Operator: A. AUTRET #53546  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	07:45:45
Diagnostic Test	Pass	07:46:20
Air Blank	0.000	07:46:47
0.100 Cal Check	0.097	07:47:08
Air Blank	0.000	07:47:38
Subject Test	0.000	07:47:56
Air Blank	0.000	07:48:24
Wait	ABT*	
Air Blank	0.000	07:49:01

\*Sequence Aborted

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001163  
Core Version: 8105.48  
05/29/2015 07:49:27

3068  
JL  
TS4

---

Standard Lot#: 11514100A2  
Last Changed By: T. GALLEGOS #28847

---

QAS: A. AUTRET #53546  
TUCSON PD

---

Operator: A. AUTRET #53546  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	07:50:21
Diagnostic Test	Pass	07:50:56
Air Blank	0.000	07:51:23
0.100 Cal Check	RFI*	07:51:36
Air Blank	IPS**	07:51:44

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001163  
Core Version: 8105.48  
05/29/2015 07:51:50

4 of 8  
JL  
TS4

---

Standard Lot#: 11514100A2  
Last Changed By: T. GALLEGOS #28847

---

QAS: A. AUTRET #53546  
TUCSON PD

---

Operator: A. AUTRET #53546  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	07:52:45
Diagnostic Test	Pass	07:53:20
Air Blank	0.000	07:53:47
0.100 Cal Check	0.096	07:54:08
Air Blank	0.000	07:54:37
Subject Test	RFI*	07:55:04
Air Blank	0.000	07:55:39

\*RFI Detect

*Incoming text message on  
cell phone during deficient  
sample*

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001163  
Core Version: 8105.48  
05/29/2015

07:56:01

*5af8*  
*MD*  
*TS4*

---

Standard Lot#: 11514100A2  
Last Changed By: T. GALLEGOS #28847

---

QAS: A. AUTRET #53546  
TUCSON PD

---

Operator: A. AUTRET #53546  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	07:57:00
Diagnostic Test	Pass	07:57:35
Air Blank	0.000	07:58:02
0.100 Cal Check	0.096	07:58:23
Air Blank	0.000	07:58:53
Subject Test	0.037*	08:02:00
Air Blank	0.000	08:02:32
Air Blank	0.000	08:03:00
Subject Test	XXX**	08:03:30
Air Blank	0.000	08:04:03

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001163  
Core Version: 8105.48  
05/29/2015 08:04:43

6 of 8  
JLD  
TSG

-----  
Standard Lot#: 11514100A2  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: A. AUTRET #53546  
TUCSON PD  
-----

Operator: A. AUTRET #53546  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	08:05:16
0.100 Cal Check	0.096	08:05:37
Air Blank	0.000	08:06:06

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

7 of 8  
TS4  
6-1-2015

QA Specialist Amy Autret 53546 AGENCY Tucson Police Department

DATE 5-29-15 TIME 0741

Intoxilyzer Serial # 80-001163 Location Lab

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 096 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE



DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001163  
Core Version: 8105.48  
05/29/2015 07:41:26

8018  
JL  
TS4

-----  
Standard Lot#: 11514100A2  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: A. AUTRET #53546  
TUCSON PD  
-----

Operator: A. AUTRET #53546  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	07:41:56
0.100 Cal Check	0.096	07:42:16
Air Blank	0.000	07:42:46

Preventative Maintenance  
Performed -- Timer Reset