



CERTIFICATION by CUSTODIAN of RECORDS
 Tucson Police Department Crime Laboratory
 1306 W. Miracle Mile, Tucson, AZ 85705

TS4
 7-2-2015

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001163 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6-25-2015 Cal 1 0.096 AC Cal 2 0.095 AC

Between the dates of 5-29-2015 and 6-25-2015 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on _____.
 The instrument was returned to service on _____.


Amy L. Autret

7/2/15

ADPS QAS Permit # 33718

SUBSCRIBED AND SWORN to before me this 2nd, day of July, 2015.

Notary Public

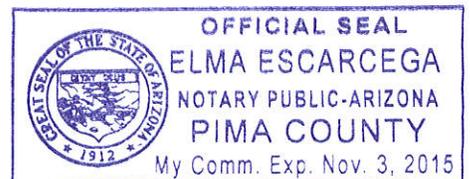

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

10F7
[Signature]
TS4
7.2.2015

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Amy Autret 53546 AGENCY Tucson Police Department

DATE 6-25-2015 TIME 1412

INTOXILYZER SERIAL # 80-001163 LOCATION Crime Laboratory

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 095 AC

Instrument is operating properly and accurately. YES NO

COMMENTS _____

SIGNATURE [Signature]

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001163
Core Version: 8105.48
06/25/2015

14:14:23

2017
Mo
TS4

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:15:22
Diagnostic Test	Pass	14:15:57
Air Blank	0.000	14:16:25
0.100 Cal Check	0.096	14:16:46
Air Blank	0.000	14:17:15
Subject Test	0.000	14:17:31
Air Blank	0.000	14:17:59
Wait	ABT*	
Air Blank	0.000	14:18:34

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001163
Core Version: 8105.48
06/25/2015 14:19:00

3 of 7
Add
TS4

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:19:53
Diagnostic Test	Pass	14:20:28
Air Blank	0.000	14:20:55
0.100 Cal Check	RFI*	14:21:11
Air Blank	IPS**	14:21:19

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001163
Core Version: 8105.48
06/25/2015

14:22:09

4 of 7
JLB

TS4

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	14:23:01
Diagnostic Test	Pass	14:23:36
Air Blank	0.000	14:24:03
0.100 Cal Check	0.096	14:24:24
Air Blank	0.000	14:24:54
Subject Test	0.203*	14:28:01
Air Blank	0.000	14:28:33
Air Blank	0.000	14:29:01
Subject Test	XXX**	14:29:22
Air Blank	0.000	14:29:56

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001163
Core Version: 8105.48
06/25/2015 14:30:12

5 of 7
JLW
B6

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:30:42
0.100 Cal Check	0.095	14:31:02
Air Blank	0.000	14:31:32

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7
TSJ
7.2.2015

QA Specialist Amy Autret 53546 AGENCY Tucson Police Department

DATE 6-25-2015 TIME 1412

Intoxilyzer Serial # 80-001163 Location Crime Laboratory

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 096 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE 

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001163
Core Version: 8105.48
06/25/2015 14:12:12

2 of 7
dd
TSh

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	14:12:41
0.100 Cal Check	0.096	14:13:02
Air Blank	0.000	14:13:31

Preventative Maintenance
Performed -- Timer Reset