



CERTIFICATION by CUSTODIAN of RECORDS
 Tucson Police Department Crime Laboratory
 1306 W. Miracle Mile, Tucson, AZ 85705

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001163 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

*TSy
8-6-2015*

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 7/23/2015 Cal 1 0.095 AC Cal 2 0.094 AC

Between the dates of 6/25/2015 and 7/23/2015 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on _____.
 The instrument was returned to service on _____.

William L Loomis

William L. Loomis


7/27/2015

ADPS QAS Permit # 33720

SUBSCRIBED AND SWORN to before me this 27th, day of July, 2015.

Elma Escarcega

Notary Public

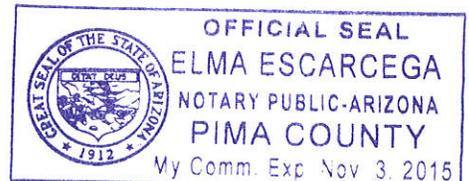

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7
Will
By
8.6.2015

QA SPECIALIST W. LOOMIS 53545 AGENCY Tucson Police Department

DATE 7/23/2015 TIME 15:21

INTOXILYZER SERIAL # 80-001163 LOCATION CRIME LAB

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 094 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Will Loomis

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001163
Core Version: 8105.48
07/23/2015 15:23:25

2 of 7
wll
TSh

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SQAP,TEST
DOB: 11/11/1987
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:24:21
Diagnostic Test	Pass	15:24:56
Air Blank	0.000	15:25:23
0.100 Cal Check	0.095	15:25:44
Air Blank	0.000	15:26:14
Subject Test	0.000	15:26:31
Air Blank	0.000	15:26:59
Wait	ABT*	
Air Blank	0.000	15:27:33

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001163
Core Version: 8105.48
07/23/2015

15:27:51

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SQAP,TEST

DOB: 11/11/1987

Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:28:40
Diagnostic Test	Pass	15:29:15
Air Blank	0.000	15:29:42
0.100 Cal Check	RFI*	15:29:57
Air Blank	IPS**	15:30:05

*RFI Detect

**Improper Sample

3 of 7
TS4
wlv

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001163
Core Version: 8105.48
07/23/2015 15:30:23

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1987
Sex: M Weight: 222

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:31:13
Diagnostic Test	Pass	15:31:48
Air Blank	0.000	15:32:15
0.100 Cal Check	0.094	15:32:36
Air Blank	0.000	15:33:06
Subject Test	0.061*	15:36:13
Air Blank	0.000	15:36:45
Air Blank	0.000	15:37:13
Subject Test	XXX**	15:37:36
Air Blank	0.000	15:38:10

*Deficient Sample

**Invalid Sample

Begin new deprivation period

4 of 7
WV
TS4

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001163
Core Version: 8105.48
07/23/2015 15:38:25

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	15:38:54
0.100 Cal Check	0.094	15:39:15
Air Blank	0.000	15:39:45

Preventative Maintenance
Performed -- Timer Reset

5 of 7
TS4
wu

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7
will
TBH 8.6.2015

QA Specialist W. LOOMIS 53545 AGENCY Tucson Police Department
DATE 7/23/2015 TIME 15:21
Intoxilyzer Serial # 80-001163 Location CRIME LAB

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 095 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Will Loomis

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001163
Core Version: 8105.48
07/23/2015 15:21:05

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: W. LOOMIS #53545
TUCSON PD

Operator: W. LOOMIS #53545
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	15:21:35
0.100 Cal Check	0.095	15:21:56
Air Blank	0.000	15:22:26

Preventative Maintenance
Performed -- Timer Reset

7 of 7 will
TS4