



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

TS6
11.25.2015

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-001163 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 11/12/15 Cal 1 0.095 AC Cal 2 0.094 AC

Between the dates of 10/15/15 and 11/12/15 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on _____.
The instrument was returned to service on _____.

[Signature]
Amy L. Autret

11/16/15


ADPS QAS Permit # 33718

SUBSCRIBED AND SWORN to before me this 15, day of November, 2015.

[Signature]

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

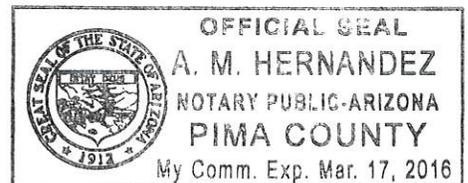

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7
TS4
11-25-2015

QA SPECIALIST Amy Autret 53546 AGENCY Tucson Police Department

DATE 11/12/15 TIME 0924

INTOXILYZER SERIAL # 80-001163 LOCATION Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS


- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 094 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

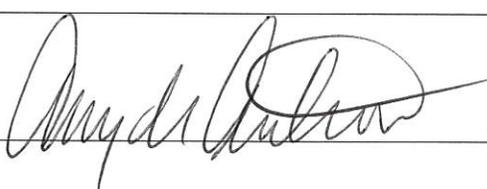SIGNATURE 

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001163
Core Version: 8105.48
11/12/2015

09:26:33

2017
TS4

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:27:27
Diagnostic Test	Pass	09:28:02
Air Blank	0.000	09:28:29
0.100 Cal Check	0.095	09:28:50
Air Blank	0.000	09:29:19
Subject Test	0.000	09:29:36
Air Blank	0.000	09:30:03
Wait	ABT*	
Air Blank	0.000	09:30:38

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001163
Core Version: 8105.48
11/12/2015 09:31:04

3067
MS
TS4

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:31:56
Diagnostic Test	Pass	09:32:31
Air Blank	0.000	09:32:58
0.100 Cal Check	RFI*	09:33:13
Air Blank	IPS**	09:33:20

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001163
Core Version: 8105.48
11/12/2015 09:33:39

4 of 7
JLB
TJ4

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:34:33
Diagnostic Test	Pass	09:35:08
Air Blank	0.000	09:35:35
0.100 Cal Check	0.096	09:35:56
Air Blank	0.000	09:36:26
Subject Test	0.166*	09:39:33
Air Blank	0.000	09:40:05
Air Blank	0.000	09:40:33
Subject Test	XXX**	09:40:55
Air Blank	0.000	09:41:28

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001163
Core Version: 8105.48
11/12/2015 09:42:02

5097
JWB
T84

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:42:31
0.100 Cal Check	0.094	09:42:52
Air Blank	0.000	09:43:21

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

057
JL
TS4
11-25-2015

QA Specialist Amy Autret 53546 AGENCY Tucson Police Department
DATE 11/12/15 TIME 0924
Intoxilyzer Serial # 80-001163 Location Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 095 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE 

DPS Form Exh G-2 (Rev. 05-01)

1/16

1/16

7 of 7
[Signature]
TS

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-001163
Core Version: 8105.48
11/12/2015 09:24:38

Standard Lot#: 11514100A2
Last Changed By: T. GALLEGOS #28847

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:25:33
0.100 Cal Check	0.095	09:25:54
Air Blank	0.000	09:26:23

Preventative Maintenance
Performed -- Timer Reset