



EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

10/17  
JLD  
754  
12.18.2015

QA SPECIALIST Amy Autret 53546 AGENCY Tucson Police Department

DATE 12-10-15 TIME 0915

INTOXILYZER SERIAL # 80-001163 LOCATION Laboratory

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS


- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 094 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE

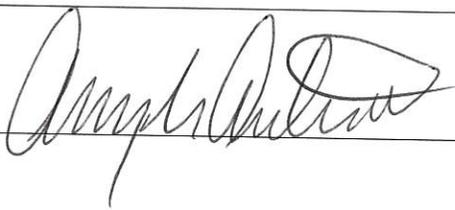

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001163  
Core Version: 8105.48  
12/10/2015

09:16:49

*2067  
Old  
TS4*

---

Standard Lot#: 11514100A2  
Last Changed By: T. GALLEGOS #28847

---

QAS: A. AUTRET #53546  
TUCSON PD

---

Operator: A. AUTRET #53546  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	09:17:41
Diagnostic Test	Pass	09:18:16
Air Blank	0.000	09:18:43
0.100 Cal Check	0.095	09:19:04
Air Blank	0.000	09:19:34
Subject Test	0.000	09:19:50
Air Blank	0.000	09:20:17
Wait	ABT*	
Air Blank	0.000	09:20:52

\*Sequence Aborted

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001163  
Core Version: 8105.48  
12/10/2015 09:21:00

3 of 7  
JLD  
TS4

---

Standard Lot#: 11514100A2  
Last Changed By: T. GALLEGOS #28847

---

QAS: A. AUTRET #53546  
TUCSON PD

---

Operator: A. AUTRET #53546  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	09:22:06
Diagnostic Test	Pass	09:22:41
Air Blank	0.000	09:23:08
0.100 Cal Check	RFI*	09:23:21
Air Blank	IPS**	09:23:30

\*RFI Detect  
\*\*Improper Sample

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001163  
Core Version: 8105.48  
12/10/2015 09:23:52

4 of 7  
old  
TS4

---

Standard Lot#: 11514100A2  
Last Changed By: T. GALLEGOS #28847

---

QAS: A. AUTRET #53546  
TUCSON PD

---

Operator: A. AUTRET #53546  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	09:24:45
Diagnostic Test	Pass	09:25:20
Air Blank	0.000	09:25:47
0.100 Cal Check	0.095	09:26:08
Air Blank	0.000	09:26:37
Subject Test	0.274*	09:29:45
Air Blank	0.000	09:30:17
Air Blank	0.000	09:30:45
Subject Test	XXX**	09:31:06
Air Blank	0.000	09:31:39

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

### Not a Successfully Completed Test Sequence

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001163  
Core Version: 8105.48  
12/10/2015 09:31:54

5067  
JL  
TS4

-----  
Standard Lot#: 11514100A2  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: A. AUTRET #53546  
TUCSON PD  
-----

Operator: A. AUTRET #53546  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:32:24
0.100 Cal Check	0.094	09:32:45
Air Blank	0.000	09:33:14

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

leaf 7  
JSH  
12-18-2015

QA Specialist Amy Autret 53546 AGENCY Tucson Police Department  
DATE 12-10-15 TIME 0915  
Intoxilyzer Serial # 80-001163 Location Laboratory

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 095 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE 

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-001163  
Core Version: 8105.48  
12/10/2015

09-15-11

*7017  
00  
TSL*

-----  
Standard Lot#: 11514100A2  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: A. AUTRET #53546  
TUCSON PD  
-----

Operator: A. AUTRET #53546  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:15:42
0.100 Cal Check	0.095	09:16:03
Air Blank	0.000	09:16:33

Preventative Maintenance  
Performed -- Timer Reset